

REGISTRATION

CANDULOR.

We have understood prosthetics since 1936.


Registration sets for full dentures,
partial and hybrid dentures, as well
as implant prosthetics.

Pocket Portfolio

The CANDULOR Pocket Portfolio offers you a
compact and clear overview of the registration
sets and their special aspects.



CRS SET 10

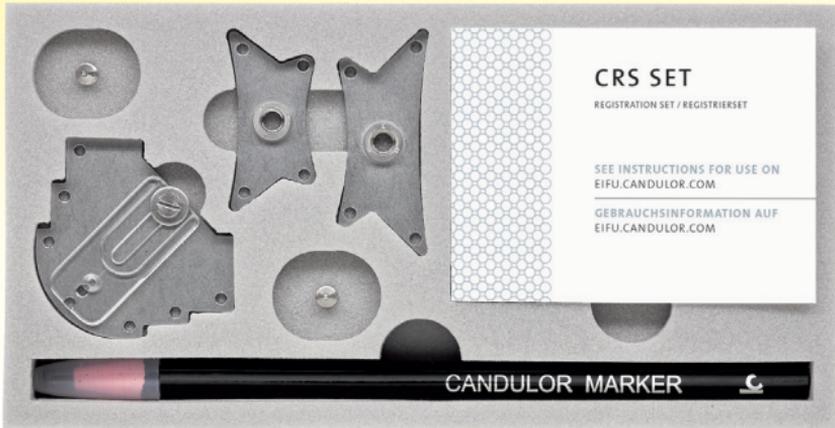

The CRS Set 10 is used for determination of the jaw relationship in complete edentulism (according to McGrane). Due to the wide plates, the set is easy to position on the edentulous jaw.

-  **Material** Stainless steel 1.4305, 1.4301 and PC
-  **Feature** Set for intraoral registration with relative bearings for full dentures
-  **Objective** Recording the position of the lower jaw in relation to the upper jaw to position the condyles centrally, laterally aligned and not displaced in the two glenoid cavities.
-  **Application** Edentulous

CRS SET 15



The CRS Set 15 is an intraoral support pin registration for the individual determination of the centric jaw position in dentulous and partially dentulous jaws and for the follow-up check of full dentures (according to McGrane).

-  **Material** Stainless steel 1.4305, 1.4301 and PC
-  **Feature** Set for intraoral registration with relative bearings for dentulous or partially dentulous jaws
-  **Objective** Determining a physiological starting point (centric)
-  **Application** Dentulous and partially dentulous


CRS SET 20



The CRS Set 20 is used for the follow-up check of full dentures and for simple determination of the jaw relationship of dentulous as well as partially dentulous jaws by means of a slow and ballistic closing movement on a frontal plateau.

-  **Material** Stainless steel 1.4305, 1.4301
-  **Feature** Set for opening and closing registration method
-  **Objective** Support of a cranioventral orientation of the condyles in the fossa (corresponds to the idealized idea of the physiological condylar position)
-  **Application** Dentulous and partially dentulous

CRS SET 10

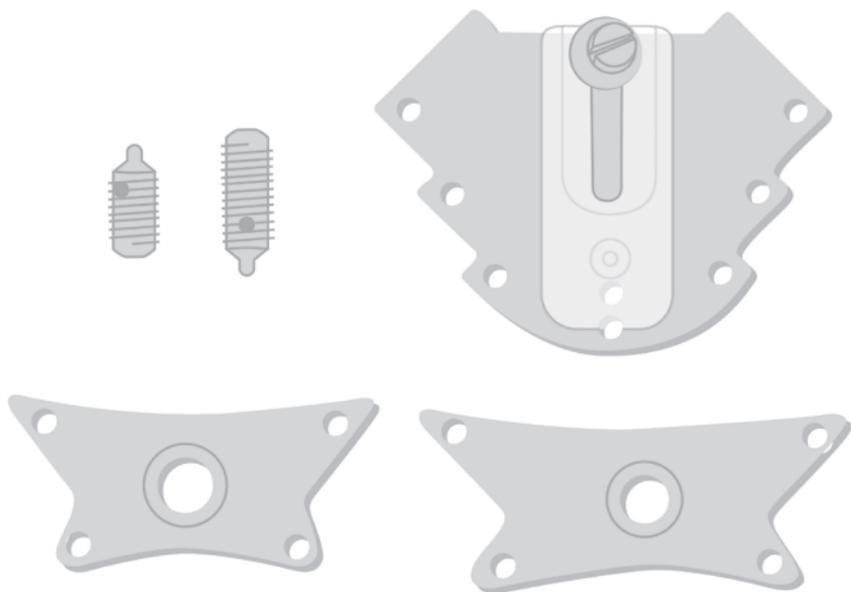

Scope of delivery

- 1 Registration tip short
- 1 Registration tip long
- 1 Registration plate
- 1 Pin supporting plate small
- 1 Pin supporting plate large
- 6 Screws and washers for the registration plate
- 6 Plexiglas plates
- 1 CANDULOR Marker

CRS SET 15



Scope of delivery

- 1 Registration plate 1
- 1 Registration plate 2
- 1 Registration plate 3
- 1 Registration plate 4
- 1 Supporting pin plate
- 1 Registration tip with sleeve short
- 1 Registration tip with sleeve long
- 6 Plexiglas roundels
- 1 CANDULOR Marker

CRS SET 20



Scope of delivery

- 2 Registration plates small
- 2 Registration plates large
- 1 Registration tip with sleeve short
- 1 Registration tip with sleeve long
- 1 CANDULOR Marker

Wish to try?

We would be pleased
to advise you.

Candolor Dental GmbH

Am Riederngraben 6


D-78239 Rielasingen-Worblingen

T +49 (0) 7731 79783-0

F +49 (0) 7731 28917

E info@candolor.de

W candolor.com

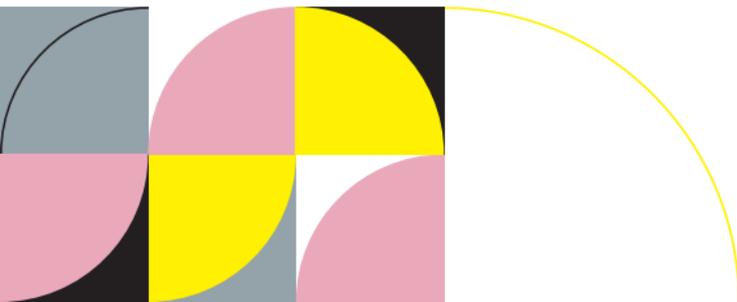

CANDULOR. HIGH END ONLY.

