



Many people say we have a think outside the box aproach, we like to think there is no box at all.

GEORGE PAPPS CEO At MPF Brush Co. our primary goal is to continue to improve our customer's performance by offering a unique range of laboratory hand tools, brushes, palettes, accessories and rotary instruments at a reasonable price, without compromising on quality.

With more than 20 years in the market we are the leaders in premium dental instruments worldwide.

Collaborating with opinion leading dental technicians from around the world has enabled us to 'think outside the box' when designing unique dental ceramic hand tools.

Thanks to our loyal customers, MPF Brush Co. has been producing dental instruments for some time now and our product range is rapidly growing. It is you that drives us to create new and innovative tools that are not only beautiful and ergonomic in design, but also functional and affordable.

We'd like to take this opportunity to thank all of our customers worldwide for their ongoing support.

1 PG 6 - 10

MPF Hybrid Series

MPF Hybrid Brush Tips (8-9)
MPF Hybrid Plus Series (10-11)

02 PG 11-33

Brushes & Brush Kits

Optimum™ Mens Brushes & Brush Kits (12)Optimum™ Ladies Brushes & Brush Kits (13)Optimum Master Brushes & Brush Kits (15)MPF Revolution Brushes & Brush Kits (16)MPF Revolution Master Brushes & Brush Kits (17)MPF Revolution Handles (19)(20)Jungo Endo Signature Brushes MPF Signature Brush (20)Stipple by Stefan Picha (20)MPF Stadard Brushes & Brush Kits (21)MPF Zirconia Stain Brushes & Brush Kit (22)MPF Synthesis Brushes & Brush Kit (23)MPF Classic Kolinsky Brushes & Brush Kit (25)(30-31)Replacement Tips

Dental Surgical

Replacement Tips

MPF Composite Stain Brushes & Brush Kit (Straight) (26)
MPF Composite Stain Brushes & Brush Kit (Angled) (26)
MPF Composite Flat Brushes & Brush Kit (Straight) (27)
MPF Composite Flat Brushes & Brush Kit (Angled) (27)
MPF Composite Application Tools Set 1 (28)
MPF Composite Application Tools Set 2 (29)

03 PG 34-42

Palettes

MPF Composite Warming Palette (35-36)
Marble Ceramic Palette (37)
Glass Ceramic Palette (38)
CAD Artistry Palette (40)
Stain/Opaque Palette (41)
Evolution Master Ceramic Palette (42)

1 PG 43-46

Carving Instruments

MPF Revolution Ceramic Sculpturing Instruments

(32-33)

05 PG 47-67

Tools & Accessories

MPF Zirco Texturizer	(48)
MPF Zir Grip	(48)
MPF Ceramo Blade	(49)
MPF Zirco Cleaner	(49)
MPF Zirco Scorer	(50)
Scribe™ Instrument by Joshua Polansky	(50)
MPF Ceramomix™ Zirconia Spatula	(51)
MPF Agate Mixing Spatula	(51)
MPF Debubblizer	(52-53)
MPF Master Klip	(54)
MPF Rainbow Drops Ceramic Colorants	(54)
MPF Brush Hydrator	(55)
MPF Brush Pool	(55)
Platinum Foil Technique by Ricardo Schäfer	(56)
MPF Drill Stands	(57)
MPF Surgical Drill Stands	(58-59)
MPF Implant Drivers	(60-61)
MPF Honeycombs & Pins	(62)
MPF Veneer & Crown Stabilizing Paste	(63)
MPF Peg Fix	(63)
MPF Waxer Instruments	(64)
MPF Crown Holders	(65)
MPF Spin Fix	(66)
Optimum™ Instrument Stand	(66)
MPF Texture Eyes	(66)
MPF All-Bright Diamond Paste	(66)
MPF Cube Stands	(67)
Liquid Dispencers	(67)
MPF Ceramic Table Mat	(67)

06 PG 68-79

Rotary Instruments


MPF ZMAX	(69)	
MPF Zirco Cera GSS	(70)	
MPF White Pearl Series	(71)	
MPF Rainbow Ceramic Polishers	(73)	
MPF Aluminium Oxide Abrasive Disc	(74)	\circ 7
MPF Polishing Brushes & Polishing Mops	(74)	O / PG 80-86
Rubber Wheels & Points: Mounted K Series	(75)	
Rubber Wheels & Points: Unmounted F Series	(76)	Dental Training School
Rubber Wheels & Points: Unmounted H Series	(77)	_
Tungsten Carbide Burs	(78)	Dental School Tools (81-84)
F.G Diamond burs	(79)	Dental School Kit (86)

Hybrid Series


MPF Hybrid Brush Tip

Patent pending (A.N.: CY202200003)


Due to the decrease of High-Quality Kolinsky, our solution was to to incorporate a Felt Nib insert within the hair bundle. This absorbent Felt Nib when wet withholds water and due to the fact that the hair bundle is wrapped around the felt nib creates a constant humid environment. This keeps the hair in a constant humid state.

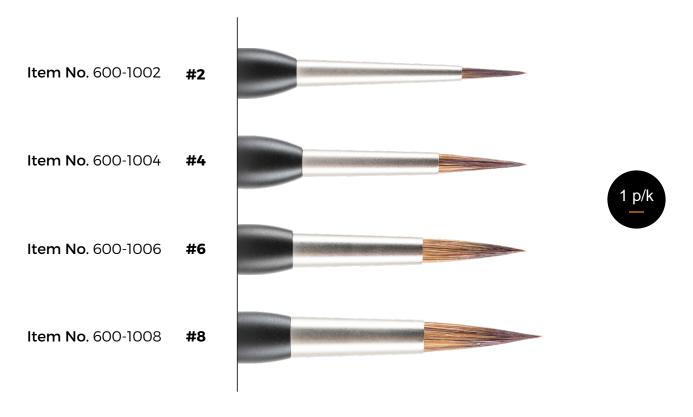

Absorbent Felt Nib

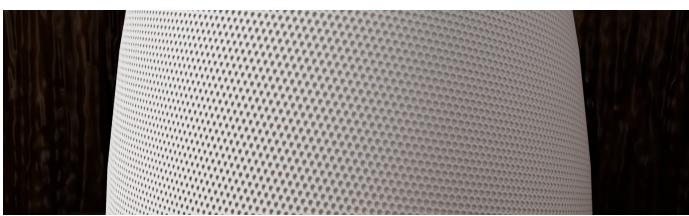

Compatible with every MPF Brush Handle.

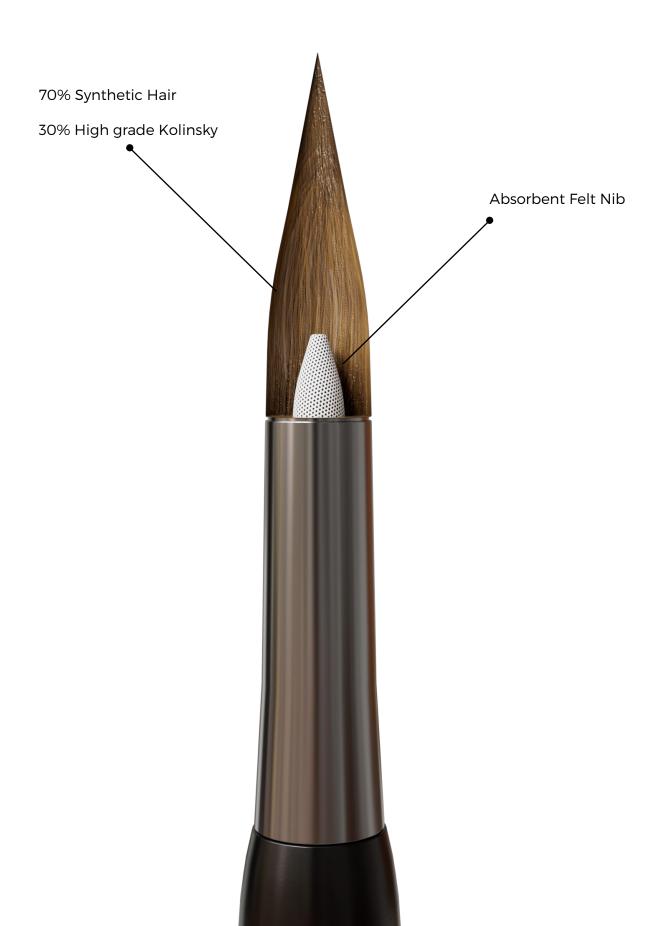

MPF Hybrid Plus Brushes & Brush Kit

Our recent Innovation in the dental brush industry, the Hybrid brush line with its absorbent internal nib implementation within the brush head, allowed us to create a much better humid environment for Synthetic hair.

The Hybrid Plus brush series are made out of 100% high grade synthetic hair that is perfect when being used for dental ceramics and staining. The absorbent material in the core of the hair bundle mimics the action of the natural hair, and delivers the exact humidity levels that are needed.

This innovative brush tip design is housed in our complete plastic handles, that also contain a few metallic balls within the ferrule of the brush tip. These balls are allowing more condensation when activated with movement.





02


Brushes & Brush Kits

Optimum™ Spring Ceramic Brushes & Kits: Mens

Spring activated tip. Perfectly balanced.Replaceable premium Kolinsky sable brush tips. Integrated thorn allows to release pinned die-segments from model base.





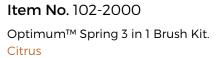



Item No. 102-3000 Optimum™ Spring 3 in 1 Brush Kit. Champagne




Item No. 102-4000 Optimum™ Spring 3 in 1 Brush Kit. Midnight Blue

Optimum™ Spring Ceramic Brushes & Kits: Ladies

Optimum™ Spring Ceramic Brushes Kits contain all 3 Brush Tip sizes (#4, #6, #8), plus the handle of your choice.





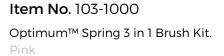



Item No. 103-2000 Optimum™ Spring 3 in 1 Brush Kit. Purple



Item No. 103-3000 Optimum™ Spring 3 in 1 Brush Kit. Midnight Blue



Optimum™ Master Brushes & Kits

All the tools you need in one package. Replaceable premium Kolinsky sable brush tips. Integrated thorn allows to release pinned die-segments from model base.

The Kit Contains: Stain, Opaque, Contour brushes, plus a Builder Brush of your choice.









Item No. 102-5006 Optimum™ Spring Master Kit. Size 6



Item No. 102-5008 Optimum™ Spring Master Kit. Size 8

MPF Revolution Brushes & Kits

Light-weight, ergonomic aluminium handle. Replaceable premium Kolinsky sable brushtips.

Revolution Same 3 tips Brush Kit contains: 3 Rplc Tips of the same size. Available size Tips: #2, #4, #6, #8

Item No. 200-0004 #4

Item No. 200-0006 #6

Item No. 200-0008 #8









Item No. 200-3002 MPF Revolution Same 3 tips Brush Kit. Midnight Black Size 2

Item No. 200-3004 MPF Revolution Same 3 tips Brush Kit. Silver Moon Size 4

Item No. 200-3006 MPF Revolution Same 3 tips Brush Kit. Lead Grey Size 6

Item No. 200-3008 MPF Revolution Same 3 tips Brush Kit. Sea Blue Size 8

MPF Revolution Master Brushes & Kits

Revolution Master Brush Kit contains: Stain, Glaze, Contour, brushes plus a Builder Brush.

Available size Tips: #2, #4, #6, #8









Item No. 200-2002 MPF Revolution Master Brush Kit. Size 2

Item No. 200-2004 MPF Revolution Master Brush Kit. Size 4

Item No. 200-2006 MPF Revolution Master Brush Kit. Size 6


Item No. 200-2008 MPF Revolution Master Brush Kit. Size 8



MPF Revolution Handles

Light-weight, ergonomic aluminium handle.

Compatible with every MPF Replacement Brush Tip.

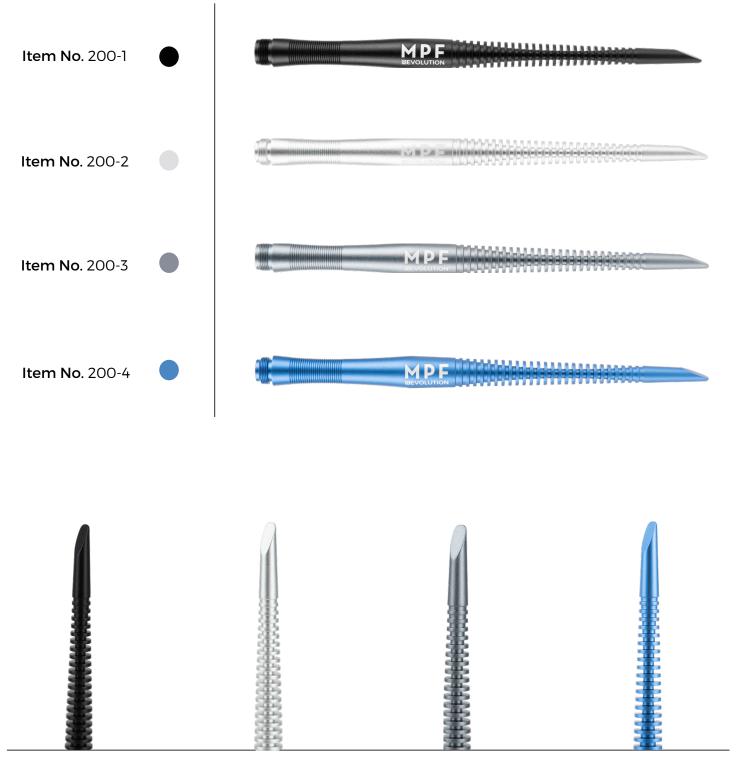

Tail Detail

Jungo Endo Signature Brushes

What makes this brush stand out from the others is its long premium kolinsky brush tip. Stay close to nature series was inspired by Jungo Endo a Japanese dental technician located in USA.



MPF Signature Brush



Stipple by Stefan Picha

A stippling brush for creating effects on the gingival surfaces.

Item No. 118-3000

MPF Standard Brushes & Kits

Innovative spring activation system. The perfect tip with a simple 'flick'. Replaceable premium Kolinsky sable brush tips.

All Standard 3 in 1 Brush Kits contain: all 3 brush tip sizes. Available sizes #4, #6, #8





Item No. 100-2000 MPF Standard 3 in 1 Brush Kit. Silver



Item No. 100-3000 MPF Standard 3 in 1 Brush Kit. Blue



Item No. 100-4000 MPF Standard 3 in 1 Brush Kit. Black

MPF Zirconia Stain Brushes & Kit

A revolutionary solution for your stain needs. Metal and rust free. Plastic handle and ferrule. Synthetic hair. Ideal for water based liquids and acids. Non-replaceable brush tips.





Item No. 105-1000 MPF Zirconia Stain 3 in 1 Brush Kit.

MPF Synthesis Brushes & Kit



The Synthesis Series Philosophy is to create a more humid environment to the synthetic hair by adapting Kolinsky into its core.


The humid properties that the kolinsky core encases are distributed through the longer synthetic hair, creating a wonderful combination at a very affordable price.





Item No. 105-0000


MPF Synthesis 4 in 1 Brush Kit

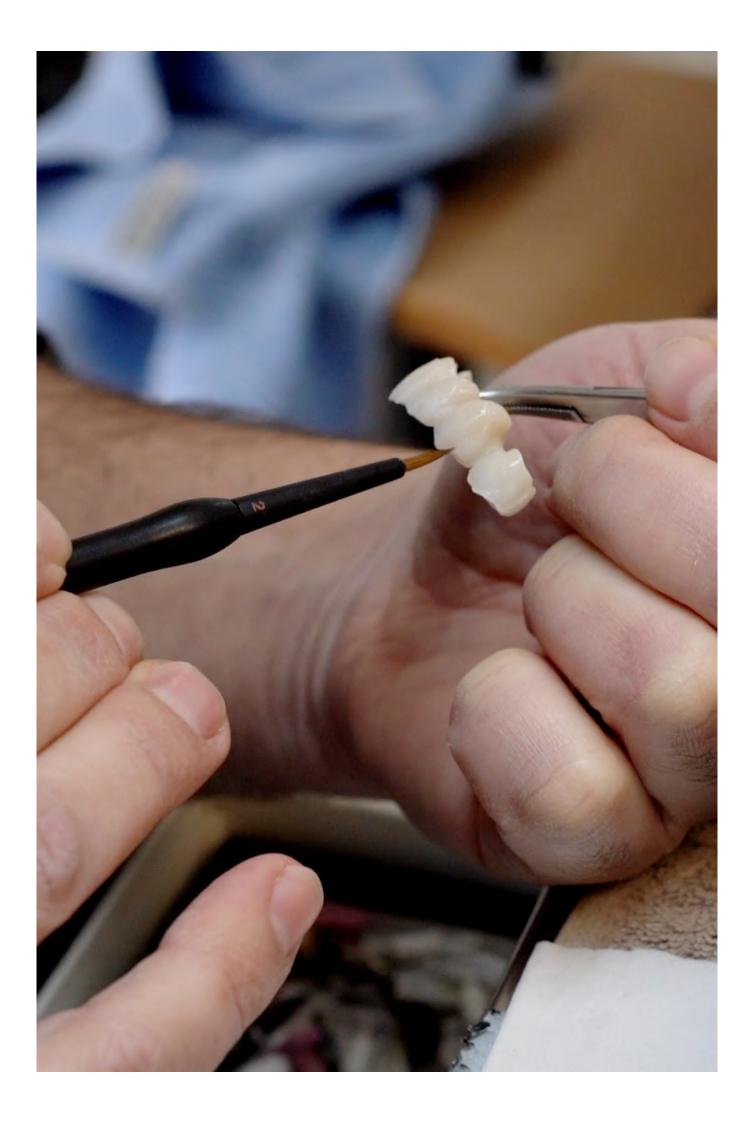

MPF Classic Kolinksy Brushes & Kit

Premium Kolinsky sable brush tips. New plastic handle and ferrule. Tip de-humidifying system (internal ball bearings.)

MPF Classic Kolinsky Brush Kit contains sizes: #2, #4, #6, #8

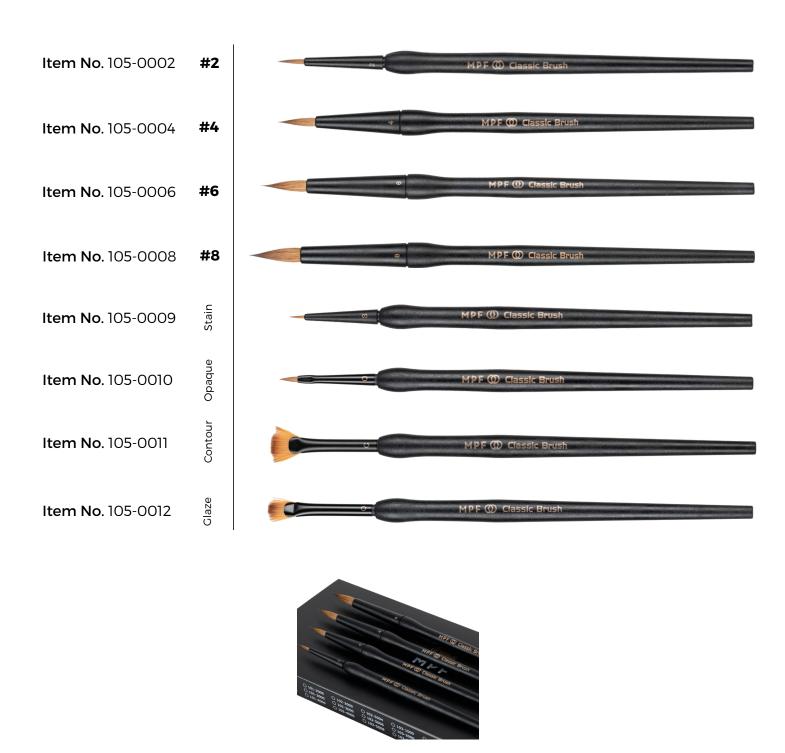


Item No. 105-0001 MPF Classic 4-1 Brush Kit

MPF Composite Stain Brushes STRAIGHT

The MPF COMPOSITE BRUSH series have been created with a special fiber optic filament that is very fine in texture but also firm to assist in application.



MPF Composite Stain Brushes ANGLED

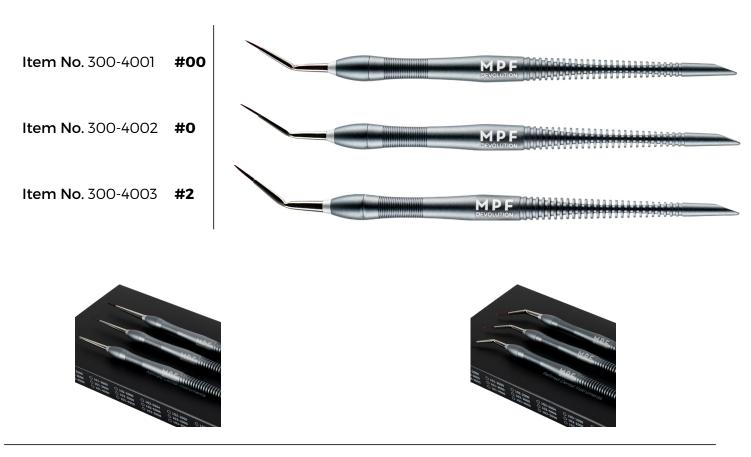

Item No. 300-3000 MPF Composite 3 in 1 Stain Brush Kit. Straight

Item No. 300-4000 MPF Composite 3 in 1 Stain Brush Kit. Angled

MPF Composite Flat Brushes STRAIGHT

Item No. 300-5001

Item No. 300-5002 #6

Item No. 300-5003 #7



MPF Composite Flat Brushes ANGLED

Item No. 300-6001 #4

Item No. 300-6002 #6

Item No. 300-6003 #7







MPF Composite 3 in 1 Flat Brush Kit. Straight



Item No. 300-6000

MPF Composite 3 in 1 Flat Brush Kit. Angled

MPF Composite Application Tools (Direct) Set 1

MPF Direct composite applicators are designed with ergonomics, comfort and balance in mind. The unique design of each and every tip allow for intraoral application and also for the dental technicians lab composite applications as well.

The handles are made from anodised aluminium and the tips are stainless steel that have a special BLACK DIAMOND coating.




Item No. 300-7000 MPF Composite Application Tools (Direct) Set 1

MPF Composite Application Tools (Direct) Set 2


These tools are all autoclavable and can withstand any sterlisation process.

They are NOT to be cleaned by using wire brush techniques as this can destroy the coating process. Care must be taken by using a soft brush to remove any reminiscence of composite material.

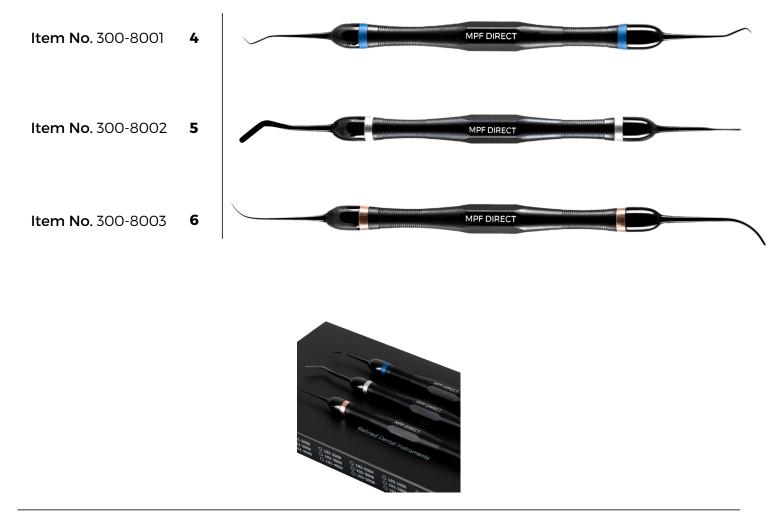

Item No. 300-8000 MPF Composite Application Tools (Direct) Set 2


Optimum™ Replacement Tips

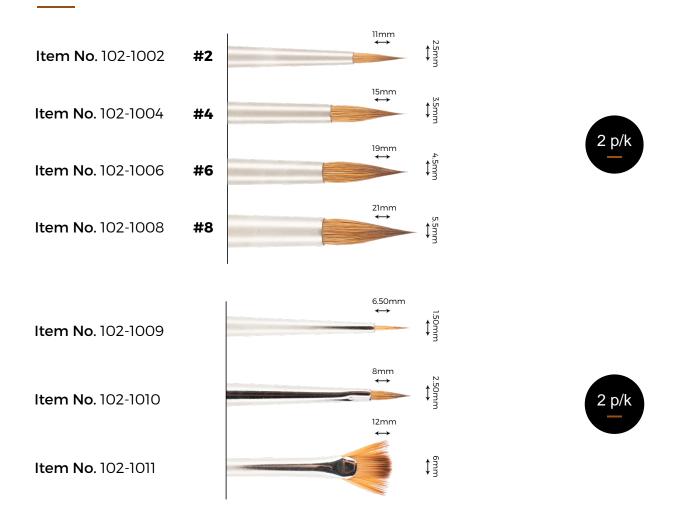


MPF Standard Replacement Tips




MPF Revolution Replacement Tips

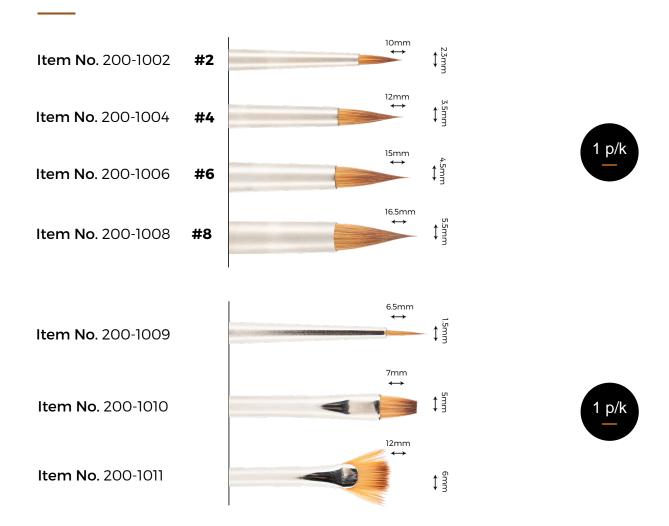

MPF Signature Replacement Tips

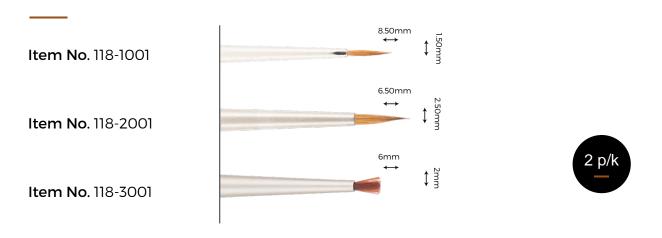

MPF Composite Stain (Straight) Replacement Tips

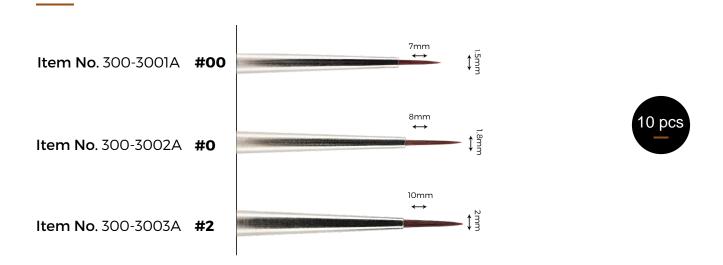


MPF Composite Stain (Angled) Replacement Tips





MPF Composite Flat (Straight) Replacement Tips

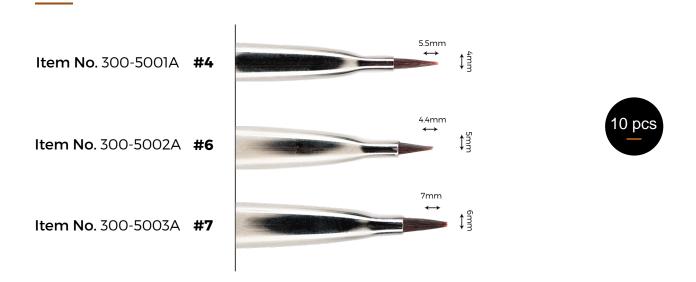

MPF Composite Flat (Angled) Replacement Tips









03

Palettes



MPF Composite Warming Palette

Introducing the first ever MPF Composite Warming Palette. This unique innovation has been designed to bring the composite to a more malleable consistency which assists the user in manipulating the composite when restoring small or long span restorations especially when the need for extensive gingival adaptation in hybrid cases is required. The palette case has been created from a polycarbonate material that is called "Makrolon", and has stable UV properties that prevents direct UV light from polymerizing the composite within the tray. Furthermore, the lid of the palette box has the ability to remain open at any height you desire, in order to prevent direct light from your laboratory surroundings.

Item No. 115-9000

Includes: Power plug cable, blade scraper.

The glass plate is warmed via an element adhered to the bottom of the glass plate which is connected via a plastic connector that is connected to the electronics on the main body of the palette. The glass plate is 6mm thick, to avoid the risk of breaking at different temperatures.

The palette features an On/Off button on the top right of the inner palette. The same button can change the temperature LCD display from Celsius to Fahrenheit by simply clicking it. When holding down the On/Off for 3 seconds the palette turns Off.

The LCD display depicts the temperature of the actual glass plate either in Celsius or Fahrenheit and can be risen or lowered via the +I - buttons. When the desired temperature has been chosen, the palette will need a few minutes in order to reach the desired temperature levels.

The palette comes with a blade scraper which is used to clear the glass plate of composite debris without scratching the glass surface. Once the composite reminisce is removed from the glass with the scraper, the glass surface can be wiped down with a clean paper towel and alcohol to clean the surface.





The palette will automatically turn off after 5 hours in case the unit has been left on by mistake.

Marble Ceramic Palette

Italian marble surface with nylon worksheet. Durable base. Capillary action keeps ceramic mix continuously wet.

Item No. 115-1000

Includes: marble tray, 2 worksheets, liquid dispenser

The unique marble used to create this palette allows controlled amounts of water to permeate and wet the nylon worksheet, ensuring a uniform, consistent and workable ceramic mix. The nylon work sheet covering the marble is kinder to Kolinsky brushes and helps keep a finer point for longer periods of time.







Item No. 115-1001 Nylon Palette Worksheets 2/pk **Item No.** 115-1002 Marble Stone

Glass Ceramic Palette

Glass work surface. Durable base. Synthetic leather strips keep ceramic mix continuously wet.


Item No. 115-0000

Includes: glass tray, spatula, 4 strips, liquid dispenser

The Glass Ceramic Palette, created for ceramists who prefer a glass working surface, uses a system of synthetic leather strips to keep ceramic mix continuously wet. Water from the base permeates up through the strips and into the ceramic stack, creating a homogenous ceramic mix.





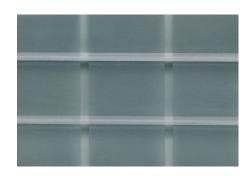

Item No. 115-0001 Palette Absorption Strips 4/pk **Item No.** 115-0002 Replacement Glass



CAD Artistry Palette

Moisture and consistency control for ceramic material. Dual working platforms. CAD-CAM Technology. This palette is the perfect solution for ceramists working with CEREC/IN-LAB and CAD-CAM technology.

Item No. 115-7000

Includes: tray, 2 strips, 1 worksheet, sponge, liquid dispenser


The CAD Artistry Palette was designed with two different options of humidification to keep the proper moisture and proper consistency of ceramic material. The two platforms on the ceramic tray allow you to choose between humidification with water or with liquid.

Choose between working on a nylon sheet or using the synthetic strips. Eleven stain compartments act as storage for your favourite stains and colours.







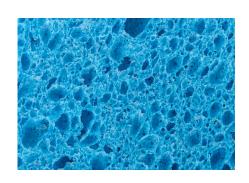

Item No. 115-7001 Nylon Palette Worksheet

Item No. 115-0001 Palette Absorption Strips 4/pk

Item No. 115-4005 Replacement Sponge

Stain/Opaque Palettes

Smart storage solution for stain and paste opaque materials. Durable plastic structure. Tight seal lid.

Item No. 115-5000 Black

Item No. 115-6000 White

The unique, angled, step-type design of the Stain & Paste Opaque Palette allows easy access to materials. Reservoirs have a very unique design, elongated towards the rear so that material flows to the back of the cup, preventing spillage. The tight seal lid prevents oxygen drying out materials. A colour chart helps keep material colours in order.



Evolution Master Ceramic Palettes

Glass work surface. Synthetic leather strips keep ceramic mix continuously wet. Concealed stain palette draw.

Item No. 115-4000

Black

Item No. 115-4100

White

Includes: glass tray, 4 strips, 1 nylon worksheet, sponge, stain brush

The Evolution Master Ceramic Palette has a sleek, contoured design. It's ergonomic design is user friendly, allowing easy access to your ceramic materials without obtrusive, robust parts inhibiting your movements to the anterior portion of the tray.

Item No. 115-4006

Replacement Glass Tray (Black)

Item No. 114-4007

Replacement Glass Tray (White)

Item No. 115-4001

Replacement Nylon Worksheet

Item No. 115-4003

Replacement Synthetic Strips 4/pk

Item No. 115-4005 Replacement Sponge





04

Carving Instruments

MPF Revolution Ceramic Sculpturing Instruments


Precision ceramic carving tools. Perfectly balanced. Replaceable tips. Ergonomic handles.

Flat Blade Item No. 200-7001 Angled Blade Item No. 200-7002 Razor Item No. 200-7003 Linqual Spoon Item No. 200-7004 12mm 5.5mm Item No. 200-1701 **Rplc Tip** 0.8mm Item No. 200-1702 **Rplc Tip** 20mm Item No. 200-1703 **Rplc Tip** Item No. 200-1704 **Rplc Tip**

Item No. 200-7000

Optimum™ Ceramic Carving Instrument

Integrated thorn allows to release pinned die-segments from model base.

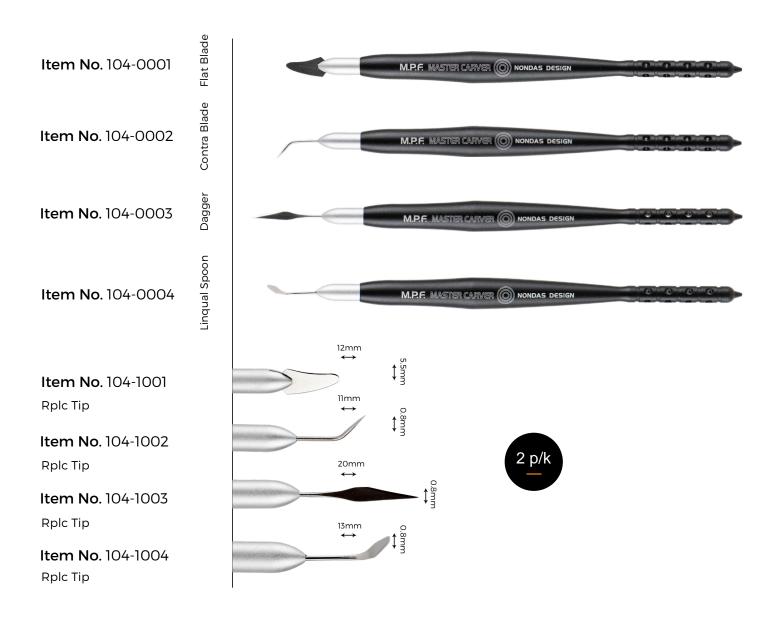



Item No. 104-3000 Optimum™ Carving Instrument Kit **Item No.** 104-2000 Optimum™ Carving Instrument Kit with Stand



M.P.F. MASTER CARVER (O) NONDAS DESIGN

m (Revers () A.M.

NASTER CARVER (Q) NONDAS DESIGN

PARINON ((C)) REVEAU PER CONTRACT

05


Tools & Accessories

MPF Zirco Texturizer

MPF Zirco Texturizer allows you to create soft suttle textures on unsintered Zirconia.

Item No. 400-1000

Before the sintering process this stainless steel wire brush can be used to create assymetric lines that look very natural after the sintering process, and also remove any fine zirconia powder that may be left on the outer surface of the restoration.

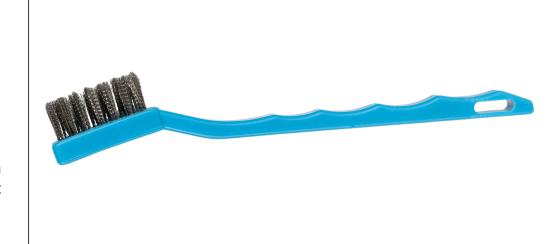

MPF Zir Grip

MPF Zir Grip Tweezers have been designed to work with all kinds of ceramic & zirconia restorations in mind.

Item No. 400-5000

These ergonomically designed tweezers boast tips made from sintered Zirconia and are shaped in such a way to assist the ceramist when positioning his crowns, veneers and small bridges onto the refractory pins before firing. These tips being from zirconia avoid leaving any dark marks that the current metallic tweezers create.



MPF Ceramo Blade

The MPF WCeramo Blade can work directly on unsintered Zirconia creating intricate surface desing.

Item No. 400-3000

Item No. 400-3001 Replacement Tip

This MPF Ceramo Blade can work directly on unsintered Zirconia creating intricate surface desing such as crack lines & deep fossa in the occlusal region, the superior sharpness allows for the ease of removing portions of zirconia to create wear facets & creative details before your final sintering.




MPF Zirco Cleaner

This newly designed brush tip has been created to allow for ease of cleaning within zirconia milled crowns, bridges & implant restoraions.

Item No. 400-4000

Item No. 400-4001 Replacement Tip

Zirconia milling creates a very fine powder after milling, and is difficult to not only see, but to remove, so this new brush tip has been designed with very stable synthetic hair in a very tightly shaped bundle, which gives excelletnt properties for removing this fine powder from your work. the fitting surface.

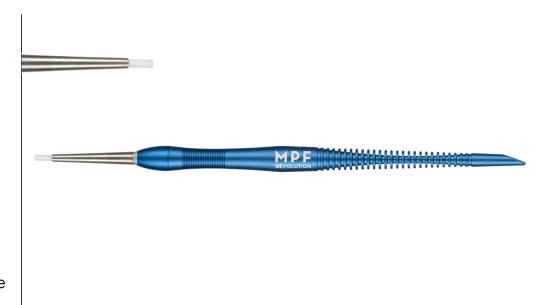

MPF Zirco Scorer

MPF Zirco Scorer allows you to create amazing surface texture on unsintered zirconia.

Item No. 400-2000

Item No. 400-2001 Replacement Tip

Due to the fact that Zirconia is a sensitive material Zirco Scorer allows us to create amazing surface texture on unsintered zirconia, pit fissure grooves and crack lines are all possible with this amazingly sharp scoring instrument and being made of Zirconia prevents any contaminations to the raw material.



Scribe™ Instrument by Joshua Polansky

Amazing surface texture and crack lines are all possible with this amazingly sharp scoring instrument.

Item No. 107-0000

Item No. 107-0001 Replacement Tip



MPF Ceramomix™ Zirconia Spatula

Ceramic mixing instrument that prevents contamination of ceramic material.

Item No. 108-1000

The spatula has been developed and combined with our spring activated technology, to relieve pressure on the zirconia blade when mixing, this has enhanced the life span of this product by 10 fold.



MPF Agate Mixing Spatula

Made from Brazilian Agate Stone, with the tail end shaped at an angle for easy extraction of ceramic powder.

Item No. 108-0000



MPF Debubblizer

The MPF Debubblizer is an advanced formula, clinically designed to reduce surface tension, thus eliminating bubble formation when pouring dental model casts from impressions and investing wax patterns. A simple mist spray into silicone, rubber base or alginate impression materials relieves surface tension that our impression materials create and this leads to a flawless replication of the impression in the thinnest and most detailed of areas. This result applies to both model replications and also wax castings.

Item No. 157-0000 300ml

Item No. 157-1000 1000ml

The MPF Debubblizer comes in a 300ml bottle and its Spray Nozzle and also 1000ml refill bottles.





Key Benefits

Reduces surface tension on wax and impressions.

Allows gypsum and investement to flow smoothly into all of your impressions.

Ensures bubble free models and smooth precision castings.

Guarantees consistent, accurate results.

MPF Master Klip

The Master Klip Tweezers are a spring activated hand tool that has stainless steel pointed holders that are activated and deactivated when the base of the tweezer is pressed.

Item No. 112-1000

These Tweezers are very useful for holding onto delicate prosthesis such as crown and veneers when staining or creating additions.





MPF Rainbow Drops Ceramic Colorants

High quality colourants for your ceramic mix .

Item No. 151-0000 (Kit)

Item No. 151-0001 (Red)

Item No. 151-0002 (Blue)

Item No. 151-0003 (Orange)

Item No. 151-0004 (Brown)

Item No. 151-0005 (Yellow)

Item No. 151-0006 (Violet)



MPF Brush Hydrator

MPF Brush Hydrator was made to keep you brushes perfectly hydrated throughout the day.

Item No. 109-0000

Pour some water in the glass reservoir, and simply mount the brush on the hydrator for about 5 minutes.

After 5 minutes of soaking, the water permeates its way into the hair bundle, giving perfect moisture level required for application.

ATTENTION: Hydrator is not a storage solution.



MPF Brush Pool

Egronomic and durable design for cleaning and storing your favourite MPF brushes. The base is made out of silicon which makes it very easy to clean and dry.

Item No. 111-0000 Brush Pool

Item No. 111-0001 Replacement Bowl



Item No. 111-0002 Replacement Sponge





Platinum Foil Technique Kit by Ricardo Schäfer

In close collaboration with one of our okey opinion leader, Ricardo Schäfer from Argentina, we developed an adapter tool that boasts a double ended tool made from aluminium that has two Zirconia adaptation tips, designed by Ricardo Schäfer



Item No. 117-0000

The Kit contains:

- Adapter Tool: these ergonomic tips have been shaped to assist in the forming of the platinum foil onto the die.
- Scissors: a pair of sharp scissors that are designed to cut the platimun foil with ease.
- Master Klip: a tool used for holding onto your dies when staining and glazing.
- Tweezers: tweezers that are designed with non invasive tips to assist in the removal of the platinum foil after the final glaze firing.



MPF Drill Stands

The MPF Drill Stand allow us to store our burs according to the task at hand. The aluminum base is quite heavy and has non-slip silicone feet for stability. Universal fit for standard Rotary Instruments & F.G Burs.





Item No. 155-0001

Item No. 155-0002





Item No. 155-0003

Item No. 155-0004



*Drills are not included

MPF Surgical Drill Stands

Our drill stand is made from Aluminum and is color anodized in three radiant exciting colors. The three level design allows for a much easier visual of the head of the drills due to the fact that they are on different levels.

On the level closest to the base you can position up to 26 friction fit burs, on the middle level you can position up to 19 friction fit burs, and on the upper level you can position up to 12 regular hand motor burs.

The Surgical drill stand is connected to a flat base that allows the top half to spin freely for quick and easy access when choosing the appropriate bur for the task at hand.

The MPF Surgical Drill stand is fully autoclave-able.


The Drill Stand should not be put into an ultrasonic with chemicals

Item No. 156-0000

Aqua

Item No. 156-0001

Citrus

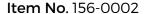

Rosegold








MPF Implant Drivers

The screwdriver bits are of the highest quality tool steel and heat treated for maximum hardness and resilience to wear and tear. The ergonomic handle is made from aluminum and color anodized and numbered in 8 brilliant colors for easy identification.

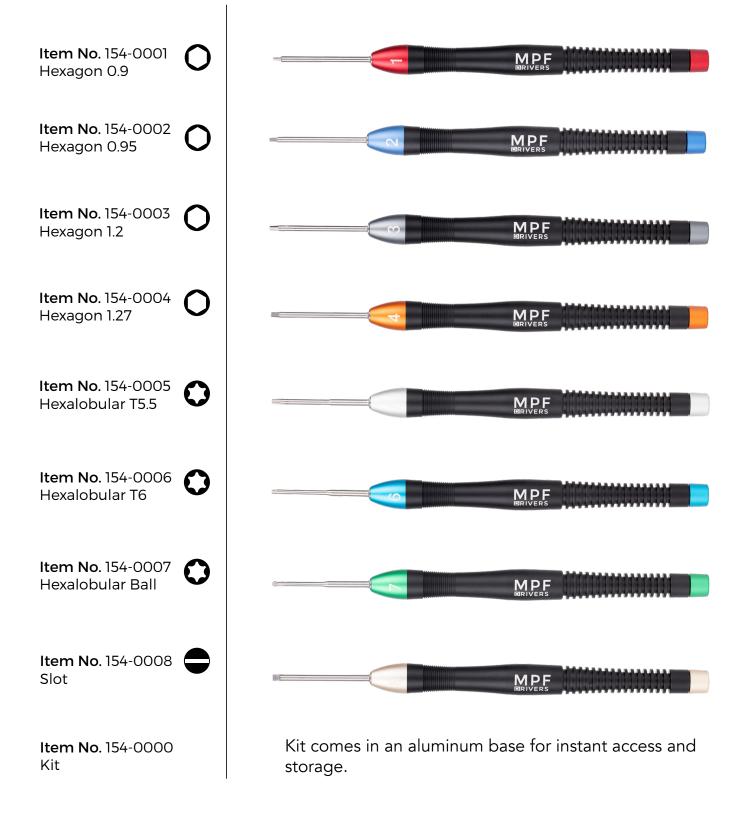






MPF IMPLANT DRIVERS

- 1. HEXAGON 0.9
- 2. HEXAGON 0.95
- 3. HEXAGON 1.2
- **4.** HEXAGON 1.27
- **5.** HEXALOBULAR T5.5
- 6. HEXALOBULAR T6
- 7. HEXALOBULAR BALL
- 8. SLOT

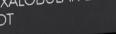

- 8. SLOT
- 7. HEXALOBULAR BALL
- 6. HEXALOBULAR T6
- 5. HEXALOBULAR T5.5
- 4. HEXAGON 1.27



MPF Honeycombs & Pins

All new 'hot zone' technology. Zirconia pins with an elevated base for crowns and implants.

Item No. 114-1000

Honeycomb Tray with 10 Ceramic Pins

Item No. 114-1001 Ceramic Pins, 10 pcs

Item No. 114-2000

Honeycomb Tray 10 Metal Pins

Item No. 114-2001

Stainless Steel Pins, 10 pcs

Item No. 114-3000

Honeycomb Tray for Implant bound CR & BR with 4 ceramic pins

Item No. 114-3001

Custom firing Ceramic Pins for Implant bound CR & CB, 4 pcs.

Item No. 114-4000

Honeycomb Tray for Veneers

40mm x 20mm

Item No. 114-4001

MPF Alumina Veneer Pins 6pcs



MPF Veneer & Crown Stabilizing Paste

A stabilizing paste used to safely adhere metal and metal free prosthesis to your model dies with superior holding properties. Burns residue free.



MPF Peg Fix

A new, revolutionary, environmentally friendly refractory stabilizing paste made from PCW fibers, suitable for firing with metal or metal free prosthesis. Can be completely removed from your prosthesis with a short blast of compressed air, leaving it totally residue free.

Item No. 150-0002 (1 tube)

Item No. 150-0003 (3 tubes)



MPF Waxer Instruments

Lightweight, aluminium handle. Replaceable stainless steel blades. Two instruments in one.

Item No. 119-0001 Spoon/Fine Dropper. Red

Item No. 119-0002 Large LeCron/Angled Spade. Blue

Item No. 119-0003 Small LeCron/Large Dropper. Purple

Item No. 119-1001 Spoon

Item No. 119-1002 Fine Dropper

Item No. 119-1003 Large LeCron

Item No. 119-1004 Angled Spade

Item No. 119-1005 Small LeCron

Item No. 119-1006 Large Dropper

Item No. 119-2000 Kit




MPF Crown Holders

Assistive tool to safely hold crowns in place while working, without chipping or marking. Special gum material is applied to the desired end (retention knurl or cut out), onto which the crown is positioned.

Item No. 152-0001 3-4mm Handle

Item No. 152-0002 2-2.5mm Handle

The MPF Crown Holder is an instrument designed to successfully hold crowns in place while working, without chipping or marking them.

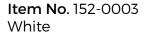

Item No. 152-0004 Black

The adhesive gum is guaranteed not to leave sticky, greasy residue on the crown or on fingers.











MPF Spin Fix



Optimum™ Instrument Stand



Item No. 106-0000

Ensures easy access and stability to crowns placed onto the tray, for ceramists to complete finishing touches with ease. Item No. 104-4000

Keeps all your brushes and tools in one place. Easy Grab-N-Go design, durable and elegant.

MPF Texture Eyes



MPF All-Bright Diamond Paste



Item No. 153-0000

An amazing copper paste that, when painted directly onto the surface of your crowns, veneers, bridge work or gypsum model, allows you to visualize and accurately replicate morphology and tooth structure.

Item No. 113-1000

Apply a small amount of paste after polishing ceramic or composite surfaces for a lustrous, high shine finish.

MPF Cube Stands





Item No. 300-1000

Ergonomic stand for storing your brushes.

Blue

Item No. 300-2000 Ergonomic stand for storing your brushes.

Black

Liquid Dispensers



Ceramic Table Mat



Item No. 115-0005 20 & 100 ml

Item No. 115-0003 **Empty Modeling Liquid**

Empty Stain Liquid Drop-

Dropper 100ml

Item No. 115-0004

per 20ml

Item No. 110-0000

06

Rotary Instruments

MPF ZMAX

Diamond impregnated discs for zirconia abrasion. Designed to generate minimal heat. RPM 12.000 - 15.000

Item No. 120-0001 Large Disc 22 x 4.5 mm

Item No. 120-0002 Small Disc 15 x 3.5 mm

Item No. 120-0003 Disc 22 x 2 mm

Item No. 120-0004 Thin Disc 22 x 0.7 mm


Item No. 120-0005 Small Barrel 4.8 x 13 mm

Item No. 120-0006 Large Barrel 6.5 x 13 mm

Item No. 120-0007 Cone 3.5 x 11 mm

Item No. 120-0008 Inverted Cone 12 x 6 mm

Item No. 120-0009 Inverted Small Cone 6 x 8 mm

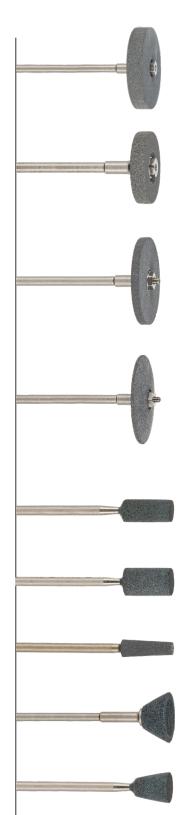

MPF Zirco Cera GSS

The GSS Series rotary instruments have extraordinary abrasive ability, they can cut through Zirconia with a minimal amount of pressure, this alleviates the chances of micro fractures within the Zirconia substructure. RPM 12.000 - 15.000

Item No. 121-0001 Large Disc 22 x 4.5 mm

Item No. 121-0002 Small Disc 15 x 3.5 mm

Item No. 121-0003 Disc 22 x 2 mm

Item No. 121-0004 Thin Disc 22 x 0.7 mm

Item No. 121-0005 Small Barrel 4.8 x 13 mm

Item No. 121-0006 Large Barrel 6.5 x 13 mm

Item No. 121-0007 Cone 3.5 x 11 mm

Item No. 121-0008 Inverted Cone 12 x 6 mm

Item No. 121-0009

Inverted Small Cone 6 x 8 mm



MPF White Pearl Series

The ultimate abrasive stones able to tackle all ceramic, zirconia and lithium di silicate materials at phenomenal speed. RPM 12.000 - 15.000

Item No. 140-0001 Large Disc 22 x 4.5 mm Item No. 140-0002 Small Disc 15 x 3.5 mm Item No. 140-0003 Disc 22 x 2 mm Item No. 140-0004 Thin Disc 22 x 0.7 mm Item No. 140-0005 Small Barrel 4.8 x 13 mm Item No. 140-0006 Large Barrel 6.5 x 13 mm Item No. 140-0007 Cone 3.5 x 11 mm Item No. 140-0008 Inverted Cone 12 x 6 mm Item No. 140-0009 Inverted Small Cone 6 x 8 mm



MPF Rainbow Ceramic Polishers

Premium diamond abrasive wheels for perfectly sealed ceramic or zirconia restoration.

Item No. 122-0001

#Coarse (Red) RPM 13000-15000

For minimal reduction on zirconia marginal regions and all types of ceramic materials.

Item No. 122-0002

#Medium (Blue) RPM 13000-15000

For minimal reduction on zirconia marginal regions and for contouring all types of ceramic materials.

Item No. 122-0003


#Fine (Grey) RPM 10000-13000

For the removal of minor scratches on all types of ceramic materials.

Item No. 122-0004

#High Shine (Yellow) RPM 8000-1000

For the final high shine procedure on all types of ceramic and zirconia materials.

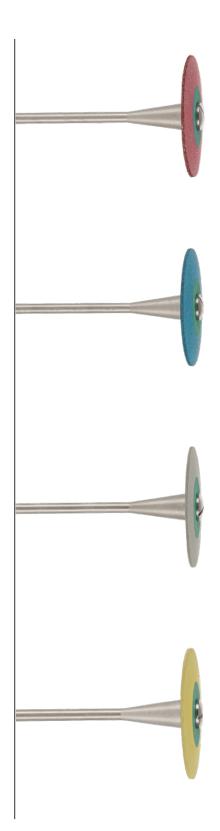

MPF Aluminium Oxide Abrasive Disc

Highly durable discs alloy restoration, without the risk of chipping or cracking the surface.

Item No. 123-0002


Large (38 x 0.6 x 1.8) 100 pcs

Item No. 123-0001

Middle (25 x 0.6 x 1.8) 100 pcs

Item No. 123-0003

Small (22 x 0.2 x 1.8) 100 pcs

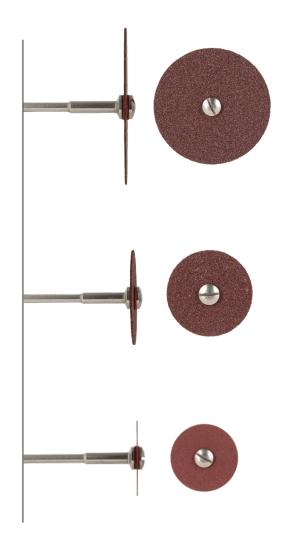

MPF Polishing Brushes & Polishing Mops

Item No. 124-0001 x100pcs

Item No. 124-0002 x50pcs

Item No. 124-0003 x25pcs

Item No. 125-0001 x100pcs

Item No. 125-0002 x50pcs

Item No. 125-0003 x25pcs



Rubber Wheels & Points. Mounted K Series.

Designed specifically to cut and polish acrylic resin base plates, denture base. Specifically for dental laboratory and clinical use.

Item No. 133-0003 (x6pcs)

Item No. 133-0001 (x25pcs)

Mounted Rubber, coarse, diam/height: 13/20

Item No. 133-0012 (x6pcs)

Item No. 133-0010 (x25pcs)

Mounted Rubber, coarse, diam/height 10/24

Item No. 133-0021 (x6pcs)

Item No. 133-0019 (x25pcs)

Mounted Rubber, coarse, diam/height 15/18

Item No. 133-0006 (x6pcs)

Item No. 133-0004 (x25pcs)

Mounted Rubber, medium, diam/height: 13/20

Item No. 133-0015 (x6pcs)

Item No. 133-0013 (x25pcs)

Mounted Rubber, medium, diam/height: 10/24

Item No. 133-0024 (x6pcs)

Item No. 133-0022 (x25pcs)

Mounted Rubber, medium, diam/height: 15/18

Item No. 133-0009 (x6pcs)

Item No. 133-0007 (x25pcs)

Mounted Rubber, fine, diam/height: 13/20

Item No. 133-0018 (x6pcs)

Item No. 133-0016 (x25pcs)

Mounted Rubber, fine, diam/height: 10/24

Item No. 133-0027 (x6pcs)

Item No. 133-0025 (x25pcs)

Mounted Rubber, fine, diam/heightL 15/18



x6pcs x25pcs R.P.M 8000



x6pcs x25pcs R.P.M 10000



x6pcs x25pcs R.P.M 8000



x6pcs x25pcs R.P.M 8000



x6pcs x25pcs R.P.M 10000



x6pcs x25pcs R.P.M 8000



x6pcs x25pcs R.P.M 8000



x6pcs x25pcs R.P.M 10000



x25pcs R.P.M 8000

Rubber Wheels & Points. Unmounted F Series.

Silicon rubber wheels and points, designed to cut and polish semi precious & precious alloys. Specifically for dental laboratory use.

Item No. 130-0004 A (x100pcs)

Item No. 130-0004 B (x30pcs*)



Rubber Wheel, coarse, diam/height: 22/3

Item No. 130-0007 A (x100pcs)
Item No. 130-0007 B (x30pcs*)



Rubber Wheel, coarse, diam/height: 22/4

Item No. 130-0001 A (x100pcs)
Item No. 130-0001 B (x30pcs*)



Rubber Wheel, coarse, diam/height: 6.5/24

Item No. 130-0005 A (x100pcs)
Item No. 130-0005 B (x30pcs*)



Rubber Wheel, medium, diam/height: 22/3

Item No. 130-0008 A (x100pcs)

Item No. 130-0008 B (x30pcs*)



Rubber Wheel, medium, diam/height: 22/4

Item No. 130-0002 A (x100pcs)

Item No. 130-0002 B (x30pcs*)



Rubber Wheel, medium, diam/height: 6.5/24

Item No. 130-0006 A (x100pcs)

Item No. 130-0006 B (x30pcs*)



Rubber Wheel, fine, diam/height: 22/3

Item No. 130-0009 A (x100pcs)

Item No. 130-0009 B (x30pcs*)



Rubber Wheel, fine, diam/height: 22/4

Item No. 130-0003 A (x100pcs)

Item No. 130-0003 B (x30pcs*)



Rubber Wheel, fine, diam/height: 6.5/24



A: x100pcs **B**: x30pcs R.P.M 10000



A: x100pcs **B**: x30pcs R.P.M 10000



A: x100pcs **B**: x30pcs R.P.M 10000



A: x100pcs **B**: x30pcs R.P.M 10000



A: x100pcs **B**: x30pcs R.P.M 10000



A: x100pcs **B**: x30pcs R.P.M 10000



A: x100pcs **B**: x30pcs R.P.M 10000



A: x100pcs **B**: x30pcs R.P.M 10000



A: x100pcs **B**: x30pcs R.P.M 10000

Rubber Wheels & Points. Unmounted H Series.

Silicon rubber wheels and points, designed to cut and polish dental ceramics. Specifically for dental laboratory use.

Item No. 132-0004 A (x100pcs)

 (\cdot)

Item No. 132-0004 **B** (x30pcs*)

Rubber Wheel, coarse, diam/height: 22/3

Item No. 132-0007 A (x100pcs)

0

Item No. 132-0007 **B** (x30pcs*)

Rubber Wheel, coarse, diam/height: 22/4

Item No. 132-0001 A (x100pcs)



Item No. 132-0001 **B** (x30pcs*)

Rubber Wheel, coarse, diam/height: 6/22

Item No. 132-0005 A (x100pcs)
Item No. 132-0005 B (x30pcs*)



Rubber Wheel, medium, diam/height: 22/3

Item No. 132-0008 A (x100pcs)
Item No. 132-0008 B (x30pcs*)

0

Rubber Wheel, medium, diam/height: 22/4

Item No. 132-0002 A (x100pcs)

Item No. 132-0002 B (x30pcs*)



Rubber Wheel, medium, diam/height: 6/22

Item No. 132-0006 A (x100pcs)



Item No. 132-0006 **B** (x30pcs*)

Rubber Wheel, fine, diam/height: 22/3

Item No. 132-0009 A (x100pcs)



Item No. 132-0009 B (x30pcs*)

Rubber Wheel, fine, diam/height: 22/4



Item No. 132-0003 A (x100pcs)
Item No. 132-0003 B (x30pcs*)

Rubber Wheel, fine, diam/height: 6/22



A: x100pcs **B**: x30pcs R.P.M 10000



A: x100pcs **B**: x30pcs R.P.M 10000



A: x100pcs **B**: x30pcs R.P.M 20000



A: x100pcs **B**: x30pcs R.P.M 10000



A: x100pcs **B**: x30pcs R.P.M 10000



A: x100pcs **B**: x30pcs R.P.M 20000



A: x100pcs **B**: x30pcs R.P.M 10000




A: x100pcs **B**: x30pcs R.P.M 10000




A: x100pcs **B**: x30pcs R.P.M 20000

Tungsten Carbide Burs


Coarse

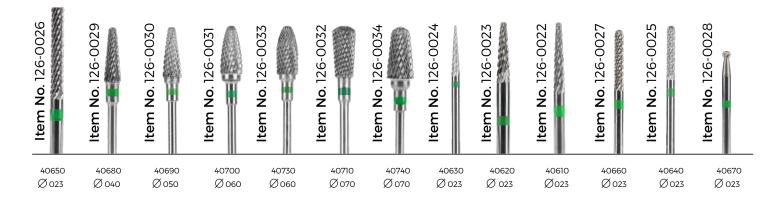


Super Coarse

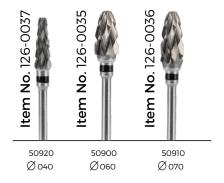

Fine

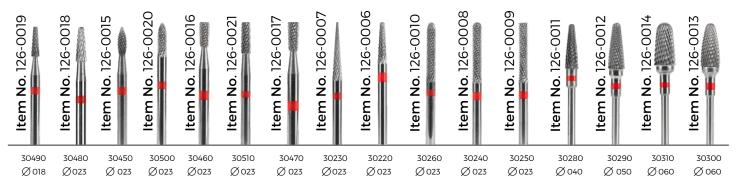



Super Fine

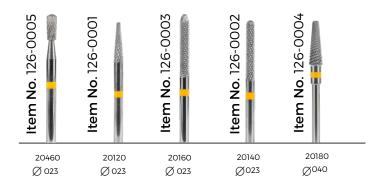

*HEAD DIAMETER: Ø

F.G Diamond Burs

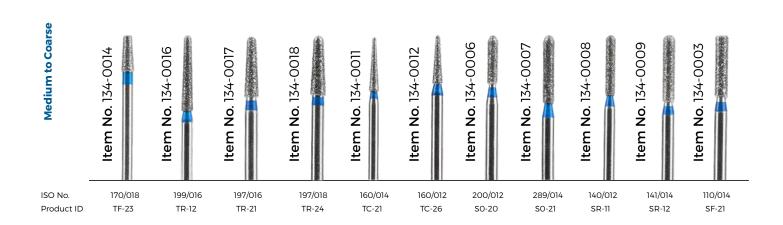

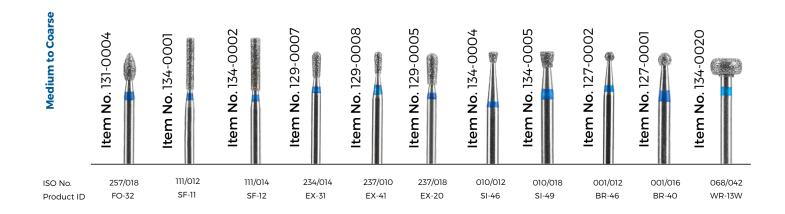





07

Dental Training School Kit

The kit comes with 27 different hand tools that will be used without a doubt in the future career of every student. We have chosen high Quality tools that will stand the test of time.

Item No. 160-0001

Large & Small wax up knive for all waxing applications, with a knife shape and also a tapered scoop shape in the one handle.

Item No. 160-0002

Item No. 160-0003

Plaster knife with a large anterior blade for cutting and a leverage shape tip at the other end of the handle, usually used when deflasking Dentures from muffles.


Item No. 160-0004

Straight edge plaster spatula for all mixing purposes.

Item No. 160-0005

Large & Small cement and acrylic mixing spatulas with a tapered tip and a flatter shaped rear tip, the handles have excellent grip due to the knurled surface of the handles.

Item No. 160-0006

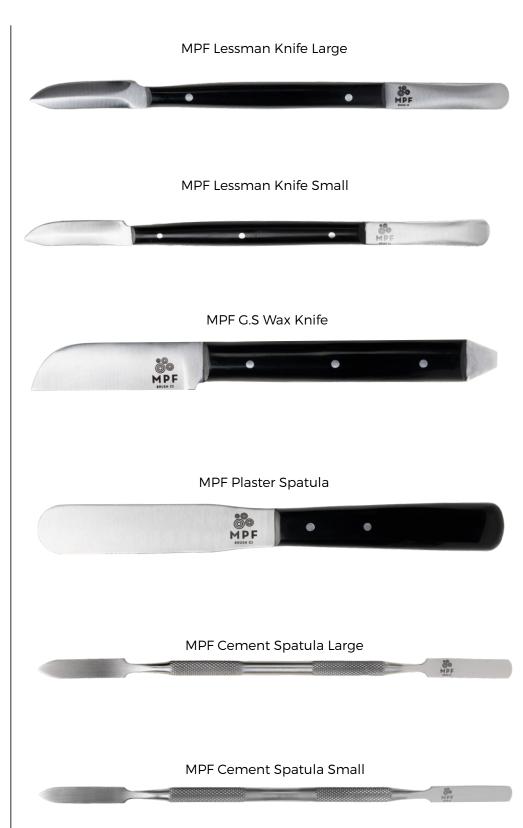

Item No. 160-0007

Microblade holder with replacement blades used for removing bubbles from plaster models and general dental model preparation.

Item No. 160-0008

Contra angle locking tweezers that can be used for holding sutures, carbon paper for bite registration purposes, general use in a Clinical or laboratory situation.

Item No. 160-0009


Number 3 Scapel handle that is used with several blade sizes for a variety of applications in a laboratory or clinical situation.

Item No. 160-0016

Straight scissors for the cutting of sutures, platinum foil and an abundance of applications in a clinical and laboratory situation.

Item No. 160-0017

Measuring calipers used for measuring the thickness of ceramic veneers and or metal framework, applied in both a clinical and laboratory situation.

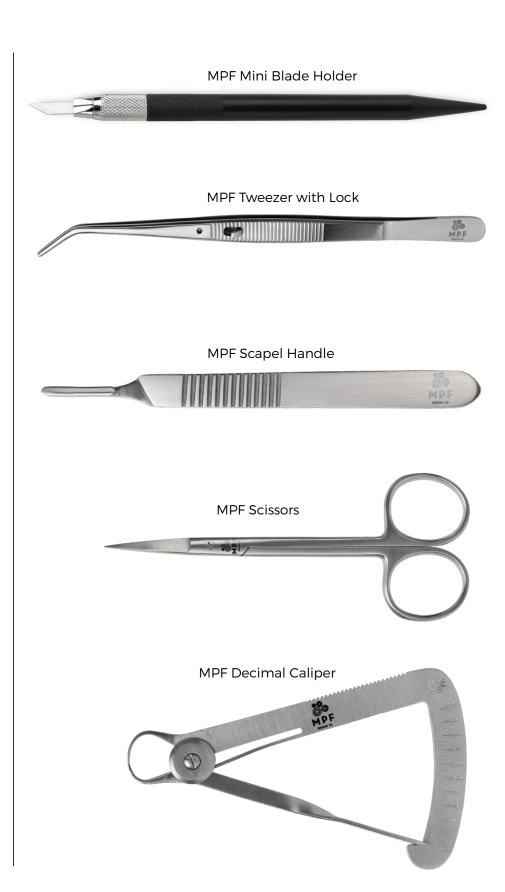

Item No. 160-0010

Concave orthodontic pliers that can be used for several applications such as creating several types of Orthodontic clasps in a Laboratory or Clinical situation.

Item No. 160-0011

Adam pliers orthodontic pliers that can be used for several applications such as creating several types of Orthodontic clasps in a Laboratory or Clinical situation.

Item No. 160-0014

Needle holder Straight can be used in surgical procedures, but also used in a dental laboratory for several applications.


Item No. 160-0015

Needle holder Curved can be used in surgical procedures, but also used in a dental laboratory for several applications.



Item No. 160-0018

Alcohol burner which is used primarily in a clinical situation when working with wax or compound materials in a surgical situation.

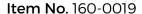

Alcohol torch can be used for the smoothening of wax surfaces in a clinical or laboratory situation.

Item No. 160-0020

Silicon Plaster or Alginate mixing bowl which can be applied in both a clinical and laboratory situation.

Item No. 160-0021

Silicon Bowls in small, medium and large sizes are primarily used for the mixing of powder to liquid materials which are applied in both a clinical and laboratory situation.



Alcohol Burner



Alcohol Torch



Plaster Bowl

Silicone Bowls







Item No. 160-0000

Dental Training School Kit *contains 27 tools in total.

Also included:

Combination Drill Kit

Ceramic Mixing Plate

MPF Classic Kolinsky Brushes #4 #6 (105-0004/105-0006)

Angle Spade/Large Lecron (119-0002)

Large Dropper/Small Lecron (119-0003)

MPF Application Tools No2 & No4 (300-7002/300-8001)





Premium Dental Instruments

MPF BRUSH Co.



MPF Brush Co.

email: info@mpfbrush.com

web: www.mpfbrush.com

tel: +35724623303