

## SILENT Extraction product navigator


Find the product that fits your needs best!

#### You haven't found the perfect Extraction yet?

Let's change that. Renfert offers suitable solutions for your requirements with its wide range of extraction systems and many years of know-how. Use this interactive guide to find the SILENT extraction system that perfectly suits you and your laboratory.

#### Work cleanly. Quietly.

Some machines do not have to be in the spotlight to make work easier. They are more like silent, constant companions. When it comes down to it, extraction units in particular, prove their value when they function without you having to pay any attention to them. Quiet, yet powerful.

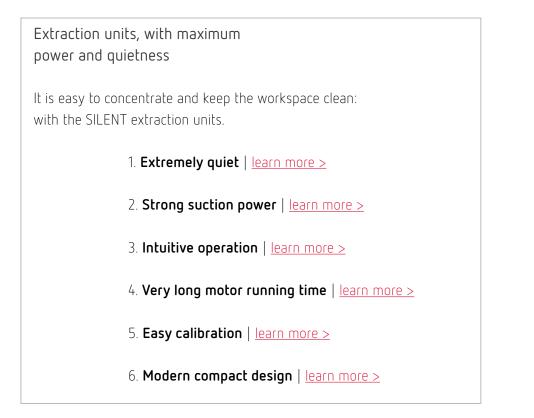

### R

#### compact



silent

#### Find your suitable product

Use our interactive product navigator to find the product you have been looking for. Or use our overview to find your extraction unit.



#### 1. Extremely quiet


#### Improve employee performance through low running noise

- Reduce employee stress-level
- Increase concentration level
- Work cleanly & quietly

One of our goals with every new extraction: as quiet as possible Some noises are more stressful than others. We are working with psychoacoustic specialists\* who will show us the unpleasant frequency ranges. This allows us to minimize the noises drastically in this range and provide the units with a pleasant audio frequency.



\*You can find more information on this topic in our white paper.

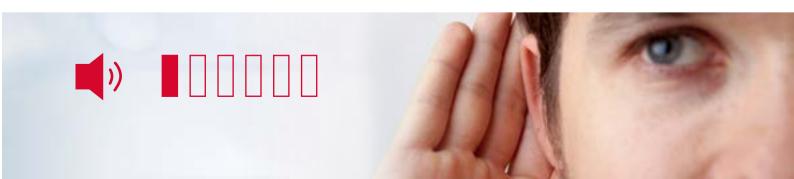

The SILENT Advantages

#### 2. Strong suction power

- Clean air for active health protection
- Dust-free workplaces
- Optimal protection of devices

#### Combine cleanliness and quietness

A workplace where the fine dust is greatly reduced and at the same time with a pleasant, low noise level greatly facilitates working. SILENT extractions all work extremely powerful with a volume flow between 2,500 and 4,000 I/min - for effective dust reduction and clean laboratory air. And always as quiet as possible.





#### 3. Intuitive operation

#### Easy start, intuitive operation

Renfert extraction units can be very easily and quickly connected (Plug & Play), operated, cleaned and maintained. The front of the unit allows convenient, direct access to all operating elements. All parameters can be very easily set.



Back

#### 4. Very long motor running time

#### Sophisticated motor technology

Two different motor technologies are available, which are both characterized by a very high service life. The collector motors provide a running time of > 1000 h. The EC motors have a long service life of > 5000 h. Fast, do-it-yourself motor change saves money and minimizes down-time.



#### Back

#### 5. Easy calibration

#### Automatic switch-on function: user-friendly

- works with ANY HANDPIECE
- works with ANY DEVICE
- works with all common CAM SYSTEMS
- COMFORTABLE & RELIABLE

#### Simply connect and start extracting

Your dust generating devices are optimally controlled together with the workstation extraction units using the easily adjustable calibration of the automatic switch-on function. The PLC interface is the intelligent bi-directional communication with your CAM unit. All you need to do is connect the extraction unit to your CAM unit, and you're ready to start.



#### 6. Modern compact design

#### Select the perfect size

SILENT extraction units are available as single-workbench, two-workbench or four-workbench extraction units and also for use with a CAD-CAM system. All versions are compact and slimline. Always compact and space-saving – it fits into every lab.



Back



Product navigator



#### Find your suitable product

Use our interactive product navigator to find the product you've been looking for. Simply choose and click on the buttons that fits your needs. Your choice will guide you to the suitable product.








#### Accessories



Product navigator > Workbench




#### Find your suitable product

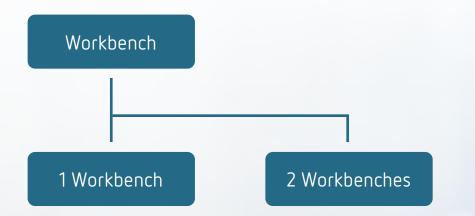

Product navigator > Workbench > 1 Workbench



#### Find your suitable product

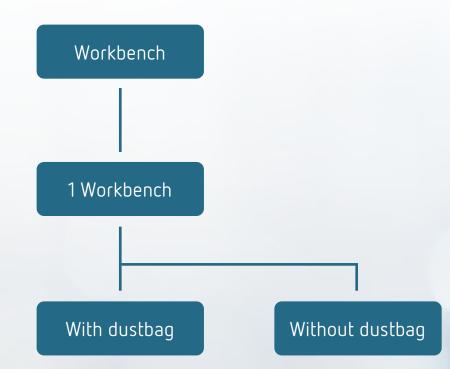












Product navigator > Workbench > 2 Workbenches



#### Find your suitable product

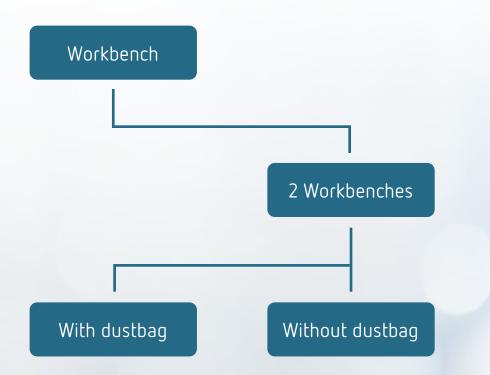







#### SILENT EC2

Two workbench extraction Powerful endurance runner with EC motor technology

Especially for strong frequented workbenches (4-8 h/day)

Visit product page



#### SILENT TC2

Two workbench extraction Performance & comfort with collector motor


For normal workloads (2-4 h/day)



Product navigator > Device




#### Find your suitable product

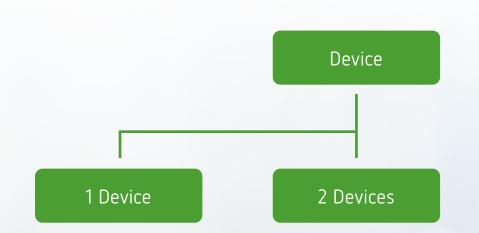

Product navigator > Device > 1 Device



#### Find your suitable product

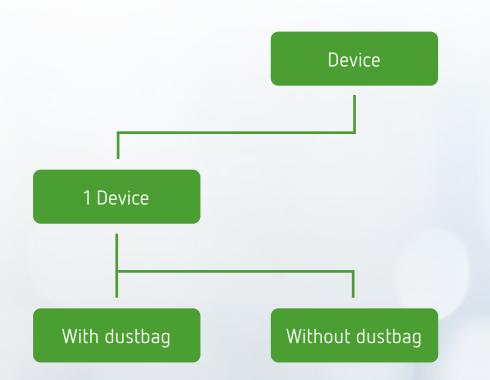



#### Vortex compact 3L

Device extraction unit Extremely powerful wet & dry extraction unit








Product navigator > Device



#### Find your suitable product

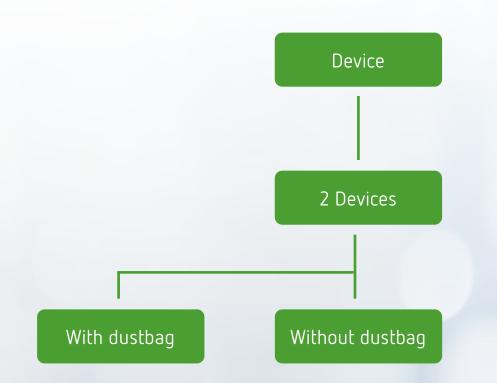



# Vortex compact 3L Device extraction unit Extremely powerful wet & dry extraction unit Visit product page Image: Ima





additionally required!

Product navigator > Device > 2 Devices > With Dustbag





additionally required!






Product navigator > CAM



#### Find your suitable product

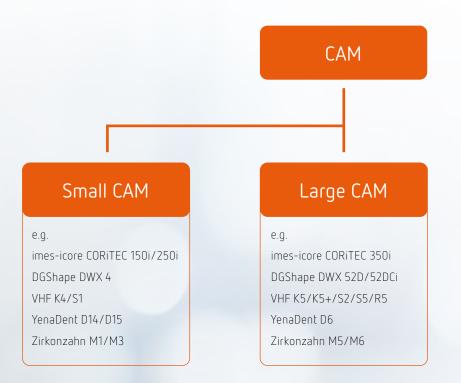





Interface cable
additionally required!





nterface cable	
additionally required!	



#### Accessoires







Suction hoses, mufflers and adapters

Visit product page









