

NEW

THE REAL



Simply **anax**, simply **THE REAL**.

THE REAL – the millable hybrid composite material


Born out of the need to find a unique product that has all the flexible properties required by frameworks without compromising on aesthetics, we came up with a Class 2a medical device that is approved for permanent use without restriction

THE REAL discs are currently made of 4-layers in 20 mm and 3-layers in 16 mm acrylic-based polymer with more than 50 % nano and micro fillers.

THE REAL – characteristics:

- 100 % biocompatible
- The monochrome or multi-layer composite blanks are characterized by their high elasticity and flexural strength.
- Gentle on abrasion and approved for final restorations due to its mechanical properties
- Multi-layer shade gradation for absolutely highly esthetic and natural restorations
- The hybrid component offers the advantage of individualization in perfection with anaxBLEND Composite Paste, Flow and anaxBLEND COLORS S
- ISO and MDR quality standards and CE-certified

THE REAL

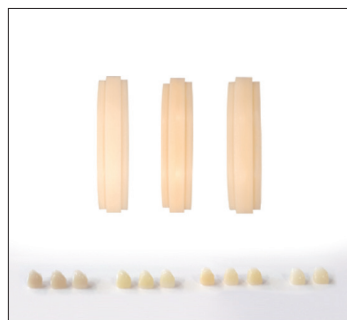

Physical properties

(based on ISO 10477:2004, Chapter5)

Flexural strength:	89 MPa (= N/mm ²)
Modulus of elasticity:	3932 MPa (= N/mm ²)
Water absorption:	9.9 ug/mm ³
Solubility:	1.5 ug/mm ³

Available in 16 mm and 20 mm height – in all standard Vita colors and on request.