

the mask range

A comprehensive range of face masks providing superior protection in a large variety of colours and patterns.



premium
plus

our type IIR face mask range

Our procedural face masks are CE marked and conform with the requirements of the Medical Devices Directive 93/42/EEC.

The Premium Plus facial protection range includes:



+ Type IIR Procedural Face Masks (Ear Loop)



+ Type IIR Procedural Face Masks (Tie-On)




+ Type IIR Procedural Face Masks with Shield (Ear Loop)

OUTER LAYER

Splash resistant non-woven fabric protective layer: effectively blocks visible substances such as droplets. The outer layer is coloured to ensure the user puts the mask on in the correct way.

MIDDLE LAYER

High density melt-blown filter layer filters potentially infective agents from the user's mouth and nose.

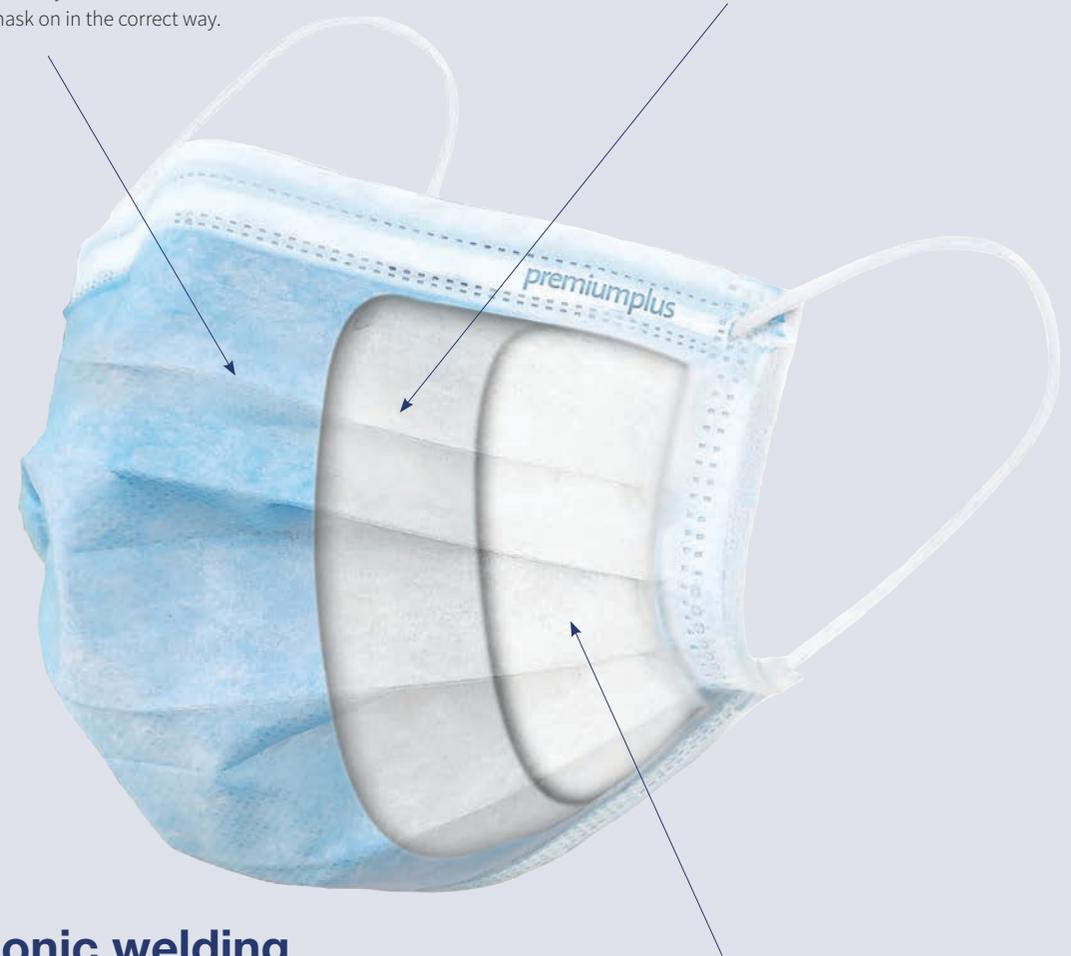

ultrasonic welding of 3 layers

3-layer filtration system is ultrasonically welded together, which provides added strength to the masks whilst ensuring a comfortable fit for the user.

INNER LAYER

Polypropylene spunbond non-woven fabric absorbs hot air exhaled from the nose and mouth, keeping the skin dry and comfortable. This ensures users can wear masks comfortably for an extended period of time.



* Approximate measurements.

4-pleat design

The 4-pleat design ensures a comfortable fit for the user by keeping space between the user and the mask, improving breathability.

BFE efficiency

Bacterial Filtration Efficiency or BFE is a measurement of the medical face mask materials ability to act as a barrier to bacterial penetration. Our face masks have a 99% BFE, which exceeds the BFE standard required under EN14683:2019 for Type IIR Masks.

latex free ear loop

Available with ties or elasticated ear loops to ensure a secure fit around the user's face.



technical data

Our face masks have been subject to thorough testing by an independent laboratory. The results of this testing is shown below:

EUROPEAN STANDARD EN14683: 2019				
	Bacterial Filtration Efficiency (%)	Differential Pressure (Pa/cm ²)	Fluid Resistance (kPa)	Microbial Cleanliness (cfu/g)
Premium Plus Type IIR	≥ 99	< 45	21.3	< 4
Type IIR	≥ 98	< 60	≥ 16.0	≤ 30
Type II	≥ 98	< 40	Testing not required.	≤ 30
Type I	≥ 95	< 40	Testing not required.	≤ 30

AMERICAN STANDARD ASTM F2100					
	Bacterial Filtration Efficiency (%)	Particulate Filtration Efficiency @0.1μ (%)	Differential Pressure (mmH ₂ O/cm ²)	Fluid Resistance (mmHg)	Flammability
Premium Plus Level 3	≥ 99	≥ 99	< 4.5	160	Class 1
Level 3	≥ 98	≥ 98	< 6.0	160	Class 1
Level 2	≥ 98	≥ 98	< 6.0	120	Class 1
Level 1	≥ 95	≥ 95	< 5.0	80	Class 1


premium plus type IIR ear loop masks

order code: 2288


Type IIR - EN14683 rated masks are high-performance face masks suitable for clinical settings to protect the patient and working environment from infective agents and to prevent cross contamination. The melt-blown polymer material used in the middle layer of the mask is a highly effective filter which prevents micro-droplets exiting the masks. Our 4-pleat design allows the wearer to comfortably expand the mask to fully cover the area from nose to chin for improved breathability. The inner layer of the masks is made of an ultra-soft non-irritating fabric, for excellent user comfort.

- + Our procedural masks adhere to European Standard EN 14683:2019 and ASTM F2100-19.
- + Available in blue, green, pink, lavender or our newest colour, black.
- + Also available with tie-on fastenings (see page 6)
- + Effective 3-layer protection designed to filter infective agents.
- + Latex free.
- + Single use only.

STANDARD COLOURS

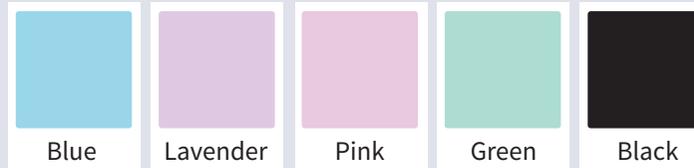

PATTERNS

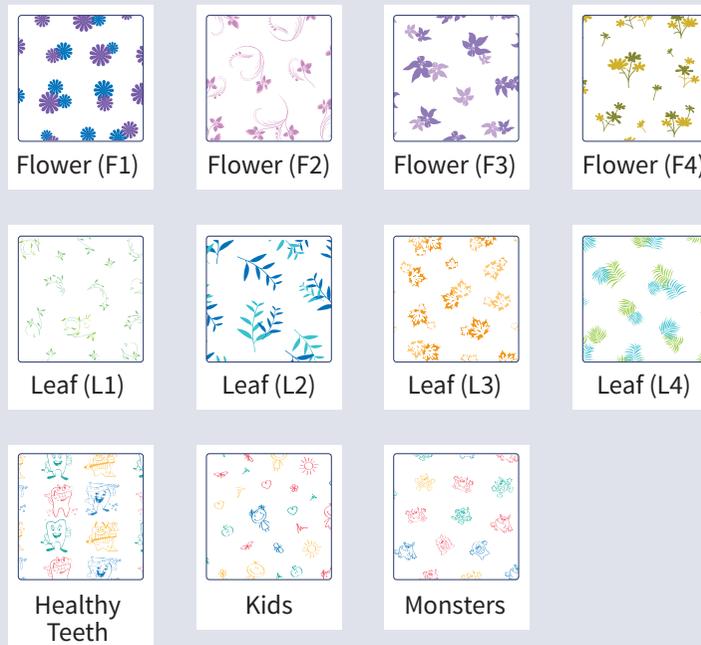

new for 2021

your favourite ear loop mask, now in black

premium plus type IIR black face masks

order code: 2288

Our best-selling ear loop mask, now available in black. This stylish mask is suitable for health professionals looking for aesthetic facial protection in clinical or other sterile environments (such as tattoo parlours or nail and beauty salons). Our new black face masks have the same high-performance qualities of our other ear loop masks, with a white inner layer to indicate to users which way the mask should be worn.



4 PLEAT DESIGN

BLACK OUTER LAYER



BLACK EAR LOOPS


- + Our procedural masks adhere to European Standard EN 14683:2019 and ASTM F2100-19.
- + Available with elastic ear loop fastenings only.
- + 4-pleat design follows the curve of the face for improved breathability and comfort.
- + Latex free.
- + Single use only.

premium plus type IIR tie-on masks

order code: 2288T

Tie-on masks are a convenient alternative to our ear loop masks. Tie-on masks are suitable for both patients and health professionals who experience discomfort or fatigue from wearing earloop masks. The 4 ultrasonically welded ties are designed to be tied horizontally behind the user's head, providing a more customised secure fitting.

- + Our procedural masks adhere to European Standard EN 14683:2019 and ASTM F2100-19.
- + Available in blue, green, pink or lavender.
- + Available with tie-on fastenings.
- + Effective 3-layer protection designed to filter infective agents.
- + Latex free.
- + Single use only.

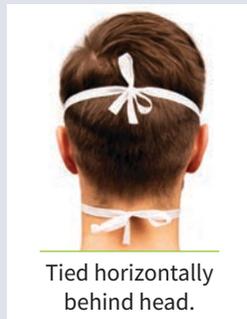

Tied horizontally behind head.

premium plus type IIR ear loop masks with shield

order code: FS2288

Our Type IIR face mask with shield provides the same high-performance protection as our normal procedural masks, with the benefit of a plastic shield which is attached to the outer surface of the mask. The shield provides added protection for the user's face and eyes against splashes and droplets.

- + Our procedure masks adhere to European Standard EN 14683:2019 and ASTM F2100-19.
- + Available in blue, green, pink or lavender.
- + Available with elastic ear loop fastenings.
- + Effective 3-layer protection designed to filter infective agents.
- + Latex free.
- + Single use only.



Shields and protects the front of the face.



STANDARD COLOURS



Blue



Lavender



Pink



Green

Order Code	Product Description	Packaging	Carton Qty
Ear Loop – Standard Colours			
2288BLUE	Face Masks Ear Loop 4 Pleat 50s Blue	50 pcs / box	20 boxes / carton
2288BLACK	Face Masks Ear Loop 4 Pleat 50s Black	50 pcs / box	20 boxes / carton
2288GREEN	Face Masks Ear Loop 4 Pleat 50s Green	50 pcs / box	20 boxes / carton
2288LAV	Face Masks Ear Loop 4 Pleat 50s Lavender	50 pcs / box	20 boxes / carton
2288PINK	Face Masks Ear Loop 4 Pleat 50s Pink	50 pcs / box	20 boxes / carton
Ear Loop - Patterns			
2288F1	Face Masks Ear Loop 4 Pleat 50s Flower F1	50 pcs / box	20 boxes / carton
2288F2	Face Masks Ear Loop 4 Pleat 50s Flower F2	50 pcs / box	20 boxes / carton
2288F3	Face Masks Ear Loop 4 Pleat 50s Flower F3	50 pcs / box	20 boxes / carton
2288F4	Face Masks Ear Loop 4 Pleat 50s Flower F3	50 pcs / box	20 boxes / carton
2288HT	Face Masks Ear Loop 4 Pleat 50s Healthy Teeth	50 pcs / box	20 boxes / carton
2288KID	Face Masks Ear Loop 4 Pleat 50s Kids	50 pcs / box	20 boxes / carton
2288L1	Face Masks Ear Loop 4 Pleat 50s Leaf L1	50 pcs / box	20 boxes / carton
2288L2	Face Masks Ear Loop 4 Pleat 50s Leaf L2	50 pcs / box	20 boxes / carton
2288L3	Face Masks Ear Loop 4 Pleat 50s Leaf L3	50 pcs / box	20 boxes / carton
2288L4	Face Masks Ear Loop 4 Pleat 50s Leaf L4	50 pcs / box	20 boxes / carton
2288MON	Face Masks Ear Loop 4 Pleat 50s Monsters	50 pcs / box	20 boxes / carton
Ear Loop with Plastic Shield			
FS2288BLUE	Face Mask Ear Loop with Shield 25s Blue	25 pcs / box	20 boxes / carton
FS2288GREEN	Face Mask Ear Loop with Shield 25s Green	25 pcs / box	20 boxes / carton
FS2288LAV	Face Mask Ear Loop with Shield 25s Lavender	25 pcs / box	20 boxes / carton
FS2288PINK	Face Mask Ear Loop with Shield 25s Pink	25 pcs / box	20 boxes / carton
Tie-On			
2288TBUE	Face Masks Tie-on 4 Pleat 50s Blue	50 pcs / box	20 boxes / carton
2288TGREEN	Face Masks Tie-on 4 Pleat 50s Green	50 pcs / box	20 boxes / carton
2299TLAV	Face Masks Tie-on 4 Pleat 50s Lavender	50 pcs / box	20 boxes / carton
2299TPINK	Face Masks Tie-on 4 Pleat 50s Pink	50 pcs / box	20 boxes / carton

order information

A Declaration of Conformity can be supplied for customers within the EU.
We are also able to supply high resolution images of our masks for promotional purposes.

product packaging

Premium Plus Type IIR Ear Loop Masks

Order code: 2288

Carton Dimensions: 53 cm x 38 cm x 22 cm (L x W x H)

Carton Weight: 4.825 kg

Cartons per EU Pallet: 28 cartons

Premium Plus Type IIR Ear Loop Masks with Shield

Order code: FS2288

Carton Dimensions: 63 cm x 38 cm x 36 cm (L x W x H)

Carton Weight: 9.370 kg

Cartons per EU Pallet: 12 cartons

Premium Plus Type IIR Tie-on Masks

Order code: 2288T

Carton Dimensions: 53 cm x 38 cm x 34 cm (L x W x H)

Carton Weight: 5.830 kg

Cartons per EU Pallet: 20 cartons

www.premiumpluspoland.com

**European Sales Office
Premium Plus Poland sp. z o.o.**

ul. Bukowska 27
62-081

Wysogotowo
Poland

+48 61 880 10 94
info@premiumpluspl.com
www.premiumpluspoland.com

**UK Sales Office
Premium Plus UK Ltd.**

Unit 2, Knighton Health Ind. Estate,
Boumemouth,
Dorset BH11 8NE

UK

+44 1202 611011
export@premiumplusuk.com
www.premiumplusuk.com

**Corporate Head Office
Premium Plus International Ltd.**

Flat 1001, Yuen Long Trading Centre,
No.33 Wang Yip Street West,
Yuen Long

Hong Kong

+852 2761 3211
mailbox@premiumplushk.com
www.premiumplushk.com

**Manufacturing Facility
Premium Plus China Ltd.**

Flat 102, No.122, Tianqiao Road
Tutang, Changping, Dongguan,
Guandong

China

+86 769 83756768
info@premiumpluschina.com
www.premiumpluschina.com

**USA Office
Premium Plus Dental Supplies Inc.**

155 Liberty Avenue
Brooklyn
New York 11212

United States of America

+1 718 366 8686
info@premiumplusdentalusa.com
www.premiumplususa.com

**Japan Office
Premium Plus Japan Co. Ltd.**

2nd Floor, 1-4-2,
Jonan, Ikedashi
Osaka, 563-0025

Japan

+81 72 735 7101
info@premiumplus.jp
www.premiumplus.jp