

MPF

BRUSH CO

PG 3-6

MPF Composite Brushes & Brush Kits

MPF Composite Stain Brushes & Brush Kit (Straight)	(3)
MPF Composite Stain Brushes & Brush Kit (Angled)	(3)
MPF Composite Flat Brushes & Brush Kit (Straight)	(4)
MPF Composite Flat Brushes & Brush Kit (Angled)	(4)
Replacement Tips	(5-6)

PG 7-10

MPF Composite Application Tools

MPF Composite Application Toos	(7-9)
MPF Composite Application Tools Kits	(10)

PG 11

Tools & Accessories

MPF Surgical Drill Stands	(11)
---------------------------	------



Video

MPF Composite Stain Brushes STRAIGHT

The MPF COMPOSITE BRUSH series have been created with a special fiber optic filament that is very fine in texture but also firm to assist in application.



MPF Composite Stain Brushes ANGLED



Item No. 300-3000

MPF Composite 3 in 1 Stain Brush Kit.
Straight



Item No. 300-4000

MPF Composite 3 in 1 Stain Brush Kit.
Angled

MPF Composite Flat Brushes STRAIGHT

Item No. 300-5001 **#4**



Item No. 300-5002 **#6**



Item No. 300-5003 **#7**



MPF Composite Flat Brushes ANGLED

Item No. 300-6001 **#4**



Item No. 300-6002 **#6**



Item No. 300-6003 **#7**



Item No. 300-5000


MPF Composite 3 in 1 Flat Brush Kit.
Straight






Item No. 300-6000

MPF Composite 3 in 1 Flat Brush Kit.
Angled

MPF Composite Stain (Straight) Replacement Tips

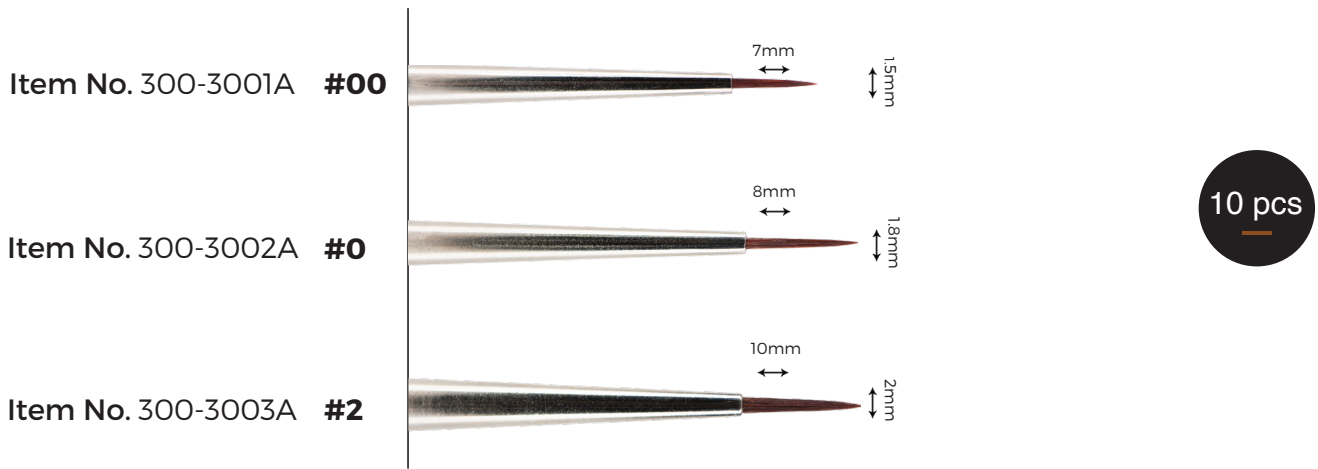


MPF Composite Stain (Angled) Replacement Tips

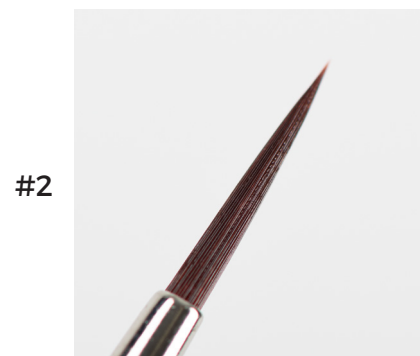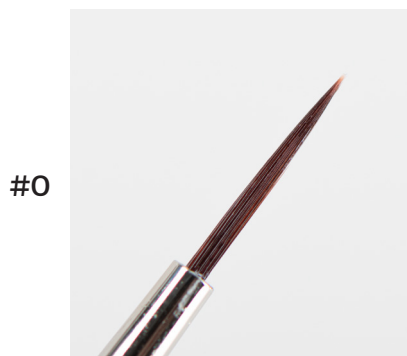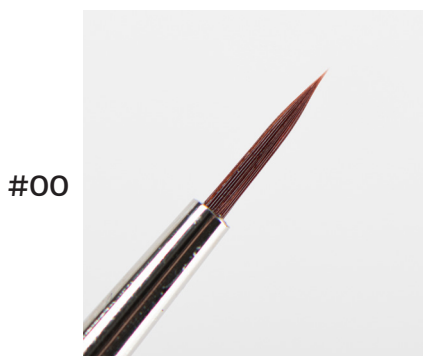

Close-up

MPF Composite Flat (Straight) Replacement Tips


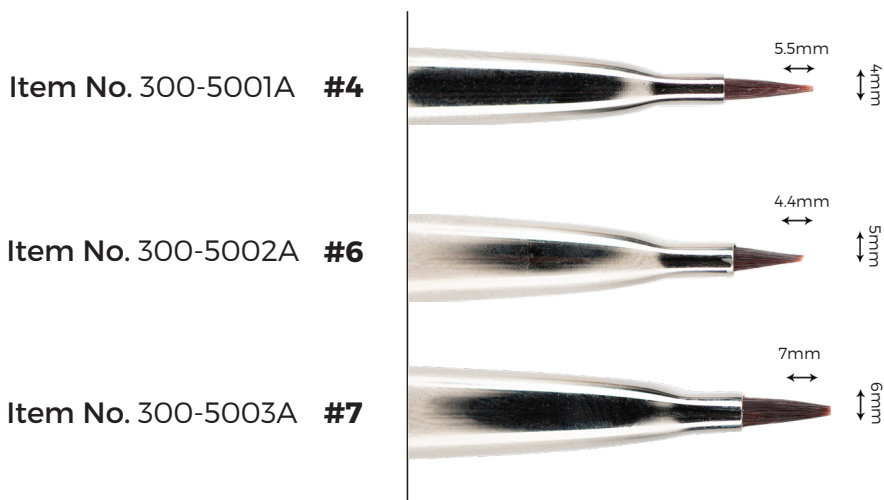


10 pcs

MPF Composite Flat (Angled) Replacement Tips



10 pcs

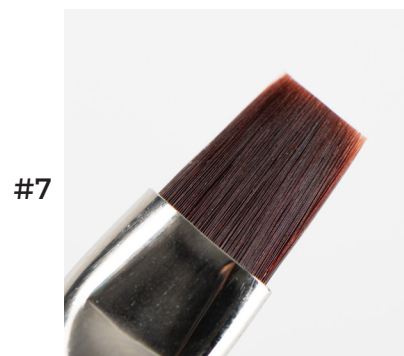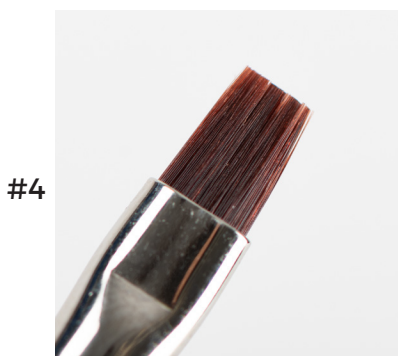

Close-up

MPF Composite Application Tools (Direct)

This Range of Composite Hand tools created by MPF have been specifically designed to enhance the users flexibility with either full mouth restorations that include gingival addition and either single crowns, inlays or onlays .

The complete range can cover all the needs that a clinician or a technician may need extraorally, from occlusal details with the finest of tips for deep fossa all the way to gingival shaping and carving and stippling effects are possible with the latest MPF composite hand tools.

Silicon tips have been added to the range to allow for ease in shaping large quantities of material such as large posterior molars and gingival addition when needed.

Assemble your tool the way it suits you:

Step 1

Chose your handle ● ● ● ——— 3 different colours for easy identification

Item No. 300-H-BL



Item No. 300-H-RG



Item No. 300-H-SL



Video

Step 2

Chose your tips



#300-7001A

kolinsky stain brush



#300-7003B

non stick silicon rounded point for gingival applications



#300-8001B

non stick silicon angled tip for the shaping of larger areas of composite



#300-8003A

synthetic stippling brush for gingival characterizations



#300-7003A

large tapered spoon shape for general application of composite in all regions



#300-7002A

flat straight tip for anterior applications in the incisal regions



#300-8001A

medium size synthetic flat brush for shaping and smoothing composite in the anterior region



#300-8002A

small tapered spoon for occlusal and incisal details



#300-7002B

angled tapered tip for interproximal applications



#300-8004B

angled spatula for mixing large qts of composite and for the shaping of gingival intended regions



#300-7001B

extra sharp deep fissure pointer for occlusal characterisations



#300-8004A

sharp deep fissure pointer for occlusal characterisations



#300-7004A

large angled occlusal groove creator



#300-7004B

long tapered deep groove and fissure creator



#300-8002B

angled deep fissure shaper for occlusal and anterior characterizations



#300-8003B

angled plunger for compressing large amounts of composite

Step 3

Assemble



MPF Composite Application Tools (Direct) Kits

Each Kit contains all 3 colors handles plus 6 specific tips.

Kit 01

#300-7000

This Kit contains:

- #300-7001A
- #300-7001B
- #300-7002A
- #300-7002B
- #300-7003A
- #300-7003B



Kit 02

#300-8000

This Kit contains:

- #300-8001A
- #300-8001B
- #300-8002A
- #300-8002B
- #300-8003A
- #300-8003B

MPF Surgical Drill Stands

Our drill stand is made from Aluminum and is color anodized in three radiant exciting colors. The three level design allows for a much easier visual of the head of the drills due to the fact that they are on different levels.

On the level closest to the base you can position up to 26 friction fit burs, on the middle level you can position up to 19 friction fit burs, and on the upper level you can position up to 12 regular hand motor burs.

The Surgical drill stand is connected to a flat base that allows the top half to spin freely for quick and easy access when choosing the appropriate bur for the task at hand.

The MPF Surgical Drill stand is fully autoclave-able.

The Drill Stand should not be put into an ultrasonic with chemicals !

Item No. 156-0000

Aqua




Item No. 156-0001

Citrus



Item No. 156-0002

Rosegold

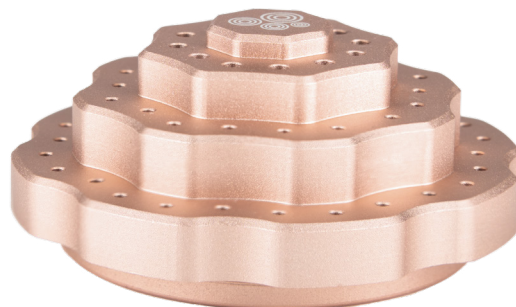

*Drills are not included

Premium Dental Instruments

MPF BRUSH Co.



MPF Brush Co. Cyprus L.T.D

web: www.mpfbrush.com

tel: +35799054381

