

MPF
BRUSH C2

“Many people say we have a think outside the box approach,
we like to think there is no box at all.”

GEORGE PAPPS
CEO

At MPF Brush Co. our primary goal is to enhance our customer's performance by offering a unique range of laboratory hand tools, brushes, palettes, accessories and rotary instruments at a reasonable price, without compromising on quality.

With over 25 years in the market we proudly stand as leaders in premium dental instruments worldwide. Our success is fuelled by collaboration with opinion-leading dental technicians from around the globe, allowing us to think outside the box and design innovative dental ceramic hand tools.

Thanks to our loyal customers, we have been producing high quality instruments for many years, and our product range continues to grow rapidly. Your feedback and support motivate us to create new and innovative tools that are not only beautiful and ergonomic in design, but also functional and affordable.

We would like to take this opportunity to extend our heartfelt thanks to all our customers worldwide for their ongoing support. Together, we will continue to shape the future of dental laboratory tools.

01

 PG 6-23

Ceramic Brushes & Brush Kits

Optimum™ Mens Brushes & Brush Kits	(7)
Optimum™ Ladies Brushes & Brush Kits	(8)
Optimum Master Brushes & Brush Kits	(9)
Optimum™ Replacement Tips	(10)
MPF Revolution Brushes & Brush Kits	(11)
MPF Revolution Master Brushes & Brush Kits	(12)
MPF Revolution Handles	(13)
MPF Revolution Replacement Tips	(14)
MPF Standard Brushes & Brush Kits	(15)
Jungo Endo Signature Brushes	(16)
MPF Standard Replacement Tips	(17)
MPF Signature Replacement Tips	(17)
MPF Classic Kolinsky Brushes & Brush Kit	(18)
MPF Xclusive Brushes & Brush Kits	(19)
MPF Hybrid Plus Brushes & Brush Kit	(20)
MPF Hybrid Brush Tips	(22)

03

 PG 33-41

Palettes

MPF Composite Warming Palette	(35)
Marble Ceramic Palette	(36)
Glass Ceramic Palette	(37)
CAD Artistry Palette	(39)
Stain/Opaque Palette	(40)
Evolution Master Ceramic Palette	(41)

05

 PG 46-50

Green State Zirconia Brushes & Tools

MPF Zirconia Brushes	(47)
MPF Zirco Texturizer	(48)
MPF Zir Grip	(48)
MPF Ceramo Blade	(49)
MPF Zirco Cleaner	(49)
MPF Zirco Scorer	(50)
MPF Ceramomix™ Zirconia Spatula	(50)

02

 PG 24-32

MPF Composite Brushes & Tools

MPF Composite Stain Brushes & Kit (Straight)	(25)
MPF Composite Stain Brushes & Kit (Angled)	(25)
MPF Composite Brushes Replacement Tips (Straight)	(26)
MPF Composite Flat Brushes & Kit (Straight)	(27)
MPF Composite Flat Brushes & Kit (Angled)	(27)
MPF Composite Brushes Replacement Tips (Angled)	(28)
MPF Composite Application Tools Direct	(29-32)

04

 PG 42-45

Paste Ceramic Characterization Kit

MPF Paste Ceramic Characterization Kits	(43)
MPF Replacement Tips	(44-45)

06 PG 51-56**Lab Materials**

MPF Veneer & Crown Stabilizing Paste	(52)
MPF Peg Fix	(52)
MPF Polishing Kit	(53)
MPF Texture Eyes	(53)
MPF All-Bright Diamond Paste	(53)
MPF Rainbow Drops Ceramic Colourants	(54)
MPF Debbublizer	(54)
MPF Honeycombs	(55)
MPF UltraSep	(56)

08 PG 70-82**Rotary Instruments**

MPF ZMAX	(71)
MPF Zirco Cera GSS	(72)
MPF White Pearl Series	(73)
MPF Rainbow Ceramic Polishers	(74)
MPF Aluminium Oxide Abrasive Disc	(75)
MPF Polishing Brushes & Polishing Mops	(75)
MPF Diamond Discs	(76)
Rubber Wheels & Points (Mounted K Series)	(77)
Rubber Wheels & Points (Unmounted F Series)	(78)
Rubber Wheels & Points (Unmounted H Series)	(79)
Tungsten Carbide Burs	(80)
F.G Diamond burs	(81)
Mandrels & Converters	(82)

07 PG 57-69**Lab Tools**

MPF Brush Hydrator	(58)
MPF Brush Pool	(58)
MPF Master Klip	(59)
Platinum Foil Technique Kit by R.Schäfer	(59)
MPF Drill Stands	(60)
MPF Surgical Drill Stands	(61)
MPF Implant Lab Screw Drivers	(62)
MPF Crown Holders	(64)
MPF Spin Fix	(65)
MPF Bench Stand	(65)
MPF Ceramic Table Mat	(66)
Liquid Dispensers	(66)
MPF Cube Stand	(66)
Scribe™ Instrument by Joshua Polansky	(67)
MPF Agate Mixing Spatula	(67)
MPF Revolution Ceramic Sculpting Instr.	(68)
MPF Waxer Tools.	(69)

09 PG 83-89**Dental Training School**

Dental School Tools	(84-87)
Dental School Kit	(89)

01

Ceramic Layering Brushes & Brush Kits

Optimum™ Spring Ceramic Brushes & Kits: Mens



Video

The Optimum™ Spring Ceramic Brushes feature a spring-activated tip that ensures perfect balance and precision. Each kit includes replaceable, premium Kolinsky sable brush tips, and an integrated thorn, on the back of the handle, designed to easily release pinned die segments from the model base.

Our kits come with all three brush tip sizes (#4, #6, #8) and your choice of handle colour, providing versatility and convenience for all your ceramic needs.



Item No. 102-2000

Optimum™ Spring 3 in 1 Brush Kit.
Citrus

Item No. 102-3000

Optimum™ Spring 3 in 1 Brush Kit.
Champagne

Item No. 102-4000

Optimum™ Spring 3 in 1 Brush Kit.
Midnight Blue

* For Replacement Tips, see page 10

Optimum™ Spring Ceramic Brushes & Kits: Ladies



Video



Item No. 103-0004 #4



Item No. 103-0006 #6



Item No. 103-0008 #8



Item No. 103-1000

Optimum™ Spring 3 in 1 Brush Kit.
Pink

Item No. 103-2000

Optimum™ Spring 3 in 1 Brush Kit.
Purple

Item No. 103-3000

Optimum™ Spring 3 in 1 Brush Kit.
Midnight Blue

* For Replacement Tips, see page 10


Optimum™ Master Brushes & Kits



Video

This kit includes replaceable premium Kolinsky sable brush tips and features an integrated thorn, on the back of the handle, for easily releasing pinned die segments from the model base.

The kit contains stain, opaque, and contour brushes, along with a building brush of your choice, ensuring you have all the essential tools for your ceramic work.

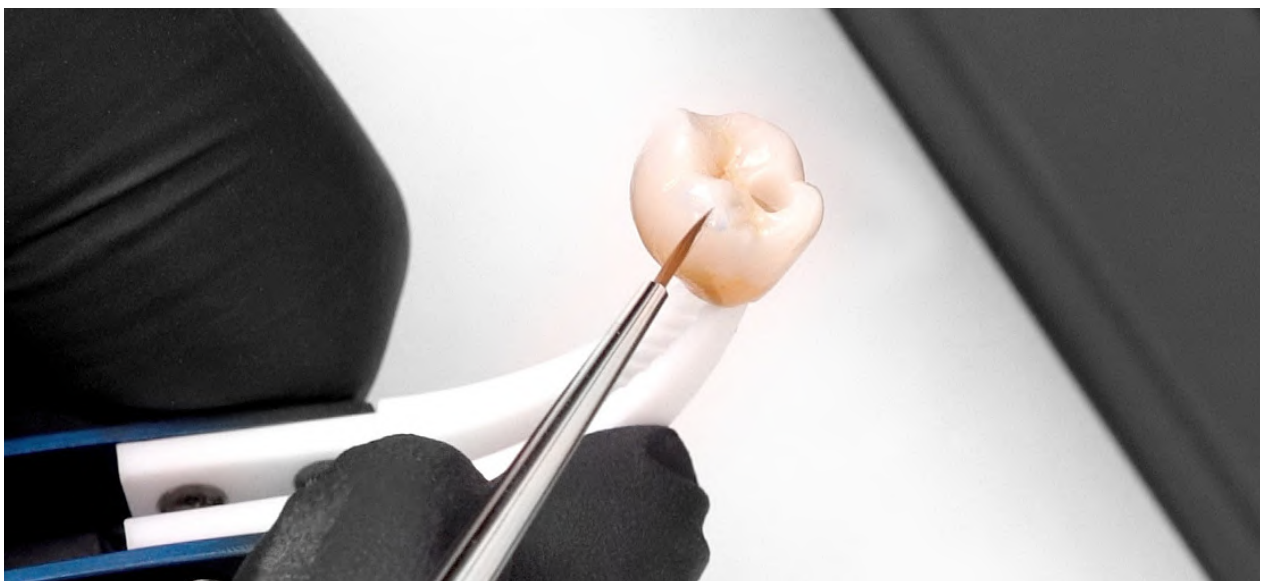


* For Replacement Tips, see page 10




Optimum™ Brushes Replacement Tips

Replaceable, premium quality tips for all your Optimum™ Brushes

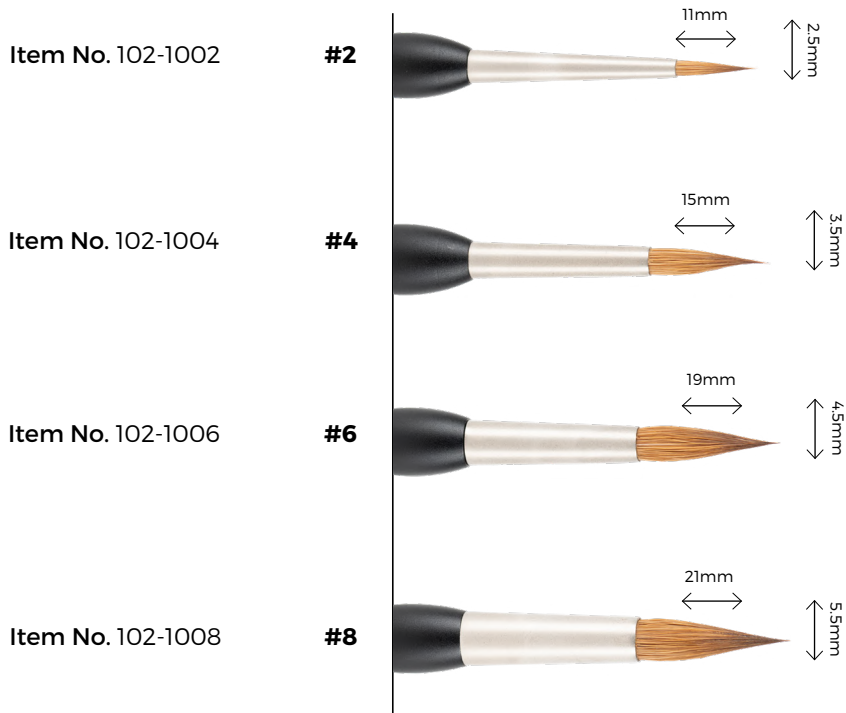

Optimum™ Master Brushes Replacement Tips

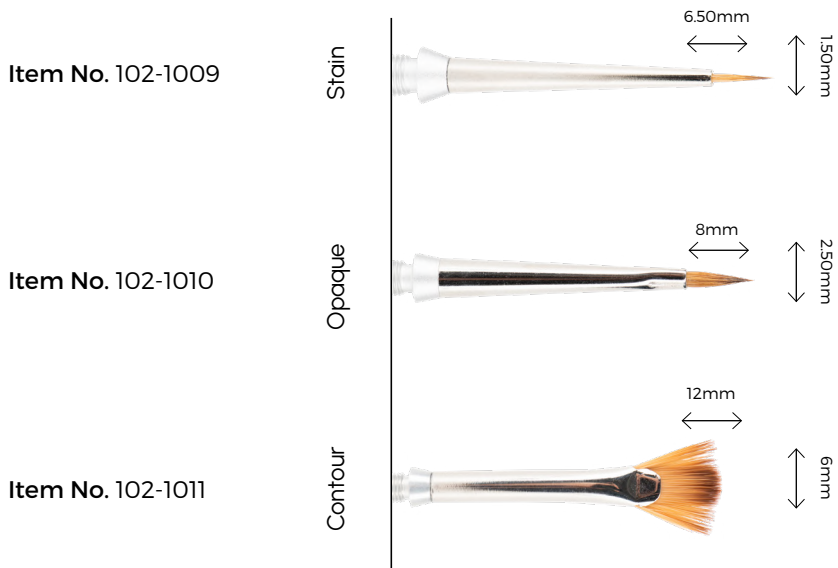

MPF Revolution Brushes & Kits



Video

Experience the perfect combination of lightweight design and ergonomic comfort with our aluminium handle. Each brush features replaceable premium Kolinsky sable tips.

The Revolution Same 3 Tips Brush Kit includes three replaceable **tips of the same size**. Available sizes are #2, #4, #6, and #8.



Item No. 200-3002

MPF Revolution Same 3 tips Brush Kit. **Midnight Black** Size 2

Item No. 200-3004

MPF Revolution Same 3 tips Brush Kit. **Silver Moon** Size 4

Item No. 200-3006

MPF Revolution Same 3 tips Brush Kit. **Lead Grey** Size 6

Item No. 200-3008

MPF Revolution Same 3 tips Brush Kit. **Sea Blue** Size 8

* For Replacement Tips, see page 14


MPF Revolution Master Brushes & Kits



Video

The Revolution Master Brush Kit includes stain, glaze, and contour brushes, along with a building brush.

Available tip sizes are #2, #4, #6, and #8.

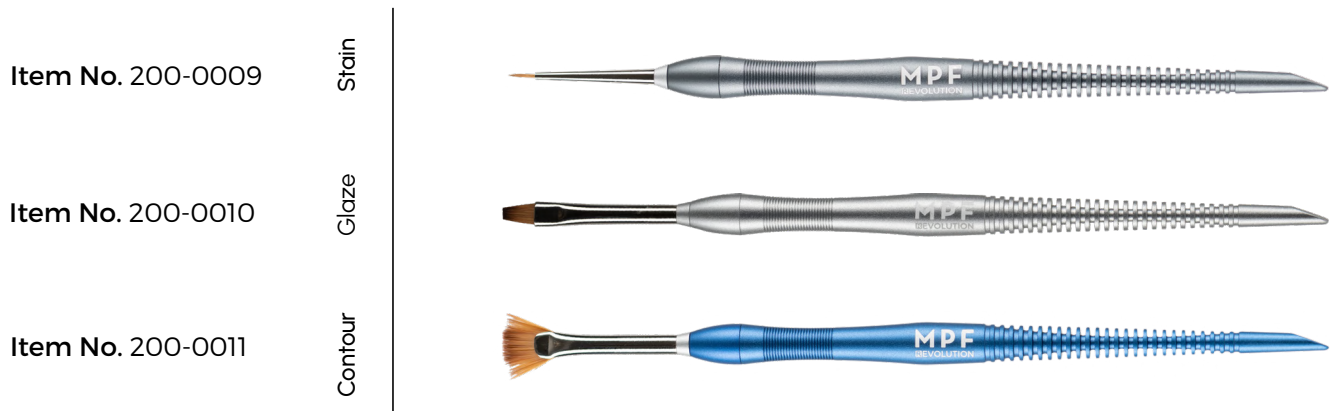

Item No. 200-2002

MPF Revolution Master Brush Kit.
Size 2

Item No. 200-2004

MPF Revolution Master Brush Kit.
Size 4

Item No. 200-2006

MPF Revolution Master Brush Kit.
Size 6

Item No. 200-2008

MPF Revolution Master Brush Kit.
Size 8

* For Replacement Tips, see page 14

MPF Revolution Handles

Light-weight, ergonomic aluminium handle.

Compatible with every MPF Replacement Builder Brush Tips.

Item No. 200-1




Item No. 200-2



Item No. 200-3



Item No. 200-4

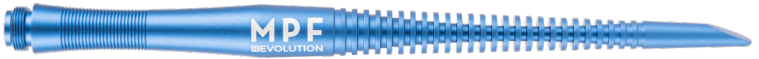


Tail Detail


MPF Revolution Replacement Tips



Replaceable, premium quality tips for all your Revolution Brushes

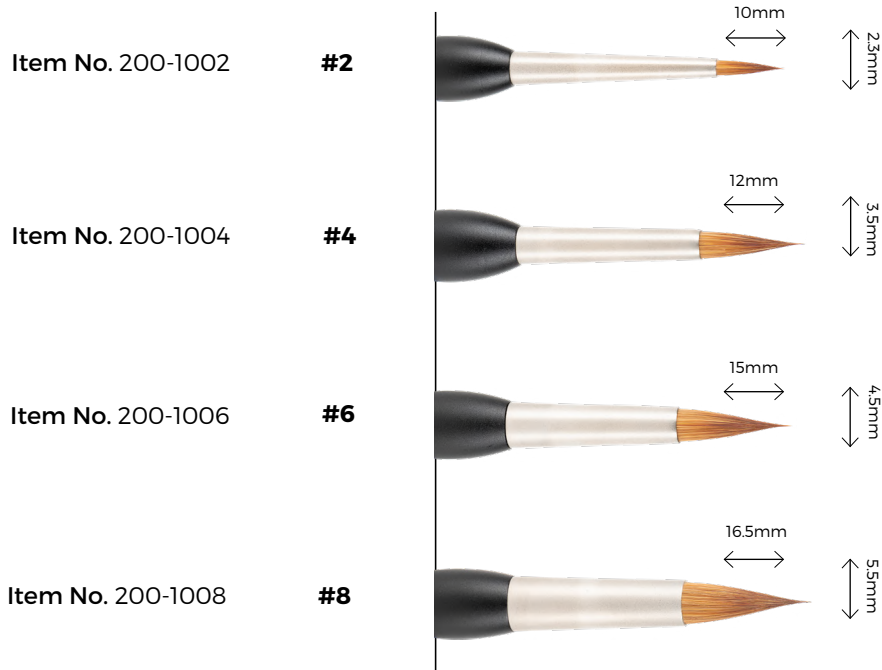

MPF Revolution Master Brushes Replacement Tips

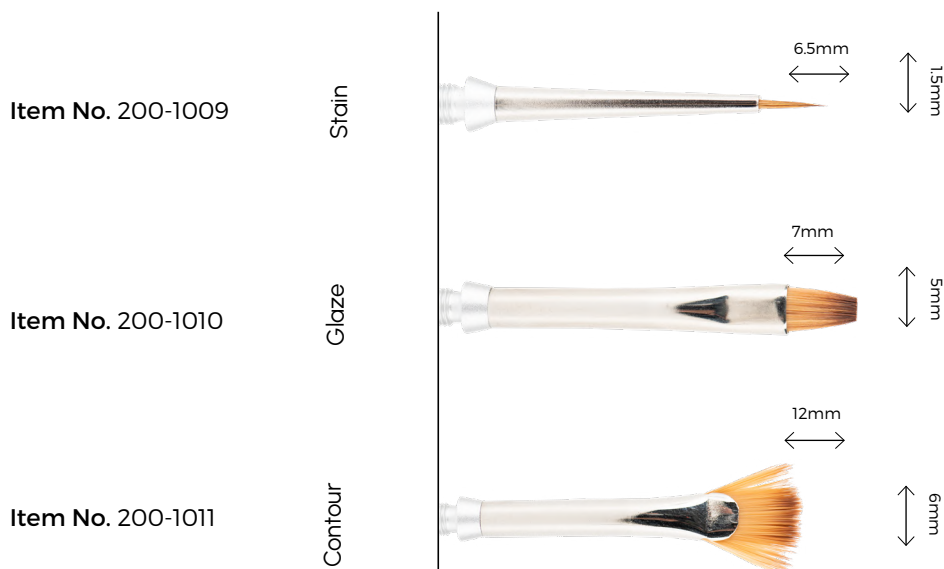





MPF Standard Brushes & Kits



Video

Experience precision and convenience with the MPF Standard Brushes, featuring an innovative spring activation system that allows for easy use with a simple flick. These brushes are designed for professionals seeking high-quality tools for their dental ceramic applications.

All Standard 3 in 1 Brush Kits come with three different brush tip sizes, providing flexibility for various tasks.

Item No. 100-0002	#2	
Item No. 100-0004	#4	
Item No. 100-0006	#6	
Item No. 100-0008	#8	



Item No. 100-2000
MPF Standard 3 in 1 Brush Kit.
Silver



Item No. 100-3000
MPF Standard 3 in 1 Brush Kit.
Blue



Item No. 100-4000
MPF Standard 3 in 1 Brush Kit.
Black

* For Replacement Tips, see page 17

Jungo Endo Signature Brushes

The Jungo Endo Signature Brush features a long premium Kolinsky sable brush tip that makes it stand out in the market.

Item No. 118-0000



Item No. 118-1000



Item No. 118-1500




* For Replacement Tips, see page 17

MPF Standard Brushes Replacement Tips



Replaceable, premium quality tips for all your Standard Brushes

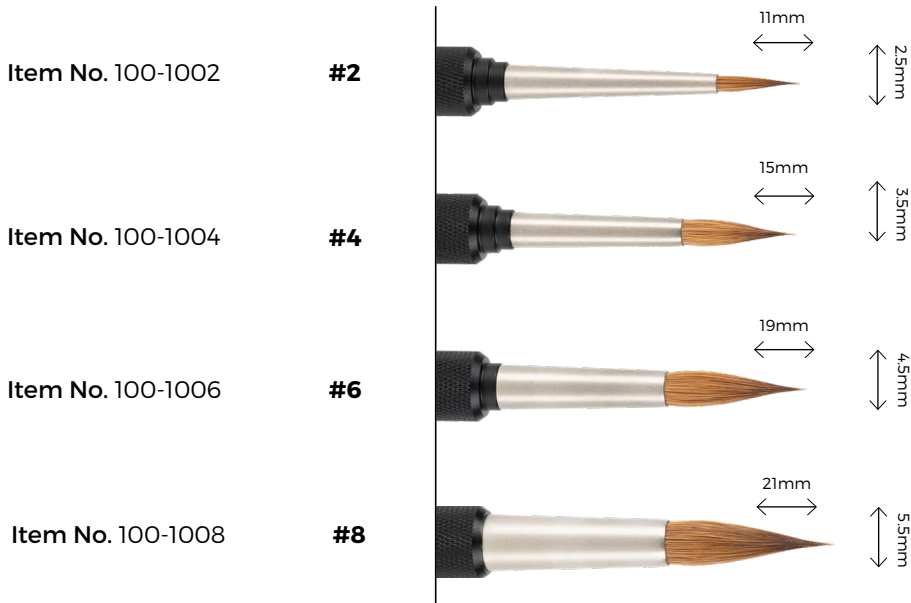

MPF Signature Brushes Replacement Tips



Replaceable, premium quality tips for all your Signature Brushes



MPF Classic Kolinsky Brushes & Kit



Video

The MPF Classic Kolinsky Brushes, expertly designed for precision and control. Featuring premium Kolinsky sable brush tips, these brushes are ideal for various dental applications.

The upgraded plastic handle ensures a lightweight, comfortable grip, while the ferrule maintains the brush's integrity and performance. Equipped with an internal ball bearing system, this innovative feature helps to maintain optimal moisture levels in the brush tips, enhancing their longevity and performance.

Item No. 105-0002	#2	
Item No. 105-0004	#4	
Item No. 105-0006	#6	
Item No. 105-0008	#8	
Item No. 105-0009	Stain	
Item No. 105-0010	Opaque	
Item No. 105-0011	Contour	
Item No. 105-0012	Glaze	



Item No. 105-0001
MPF Classic 4 in 1 Brush Kit

MPF Xclusive Kolinsky Brushes & Kit

With our extensive experience sourcing Kolinsky hair from various regions, we are able to enhance production efficiency allowing us to offer the MPF XCLUSIVE Brush at a competitive price, making it an ideal choice for markets that prioritize affordability without sacrificing performance. This brush features fine 100% natural Kolinsky tails blended at varying lengths to create a stiffer hair bundle

The lightweight, recyclable Birch wood handle is beautifully lacquered, giving each brush a unique aesthetic appeal. Its ergonomic design includes an angulated tail end, which helps effortlessly remove ceramic from bottles, while the flat end serves as a practical aid during ceramic construction.

Item No. 600-0002 **#2**



Item No. 600-0004 **#4**



Item No. 600-0006 **#6**



Item No. 600-0008 **#8**



Item No. 600-0000

MPF Xclusive 4 in 1 Brush Kit



MPF Hybrid Plus Brushes & Brush Kit

Our innovation in the dental brush industry, the Hybrid Brush Line, features an absorbent internal nib integrated into the brush head, creating an improved humid environment for synthetic hair. The Hybrid Plus brush series is crafted from 100% high-grade synthetic hair, making it ideal for use with dental ceramics and staining.

The absorbent core material mimics the properties of natural hair, providing the precise humidity levels required for optimal performance. This innovative brush tip design is complemented by durable and recycled plastic handles, which incorporate metallic balls within the ferrule. These balls enhance condensation when the brush is in motion, incorporates further functionality. The brush can be fully recycled at the end of the product cycle.


Item No. 700-0009	Stain	
Item No. 700-0002	#2	
Item No. 700-0004	#4	
Item No. 700-0006	#6	
Item No. 700-0008	#8	



Item No. 700-0000

MPF Hybrid Plus 4 in 1 Brush Kit

100% Synthetic

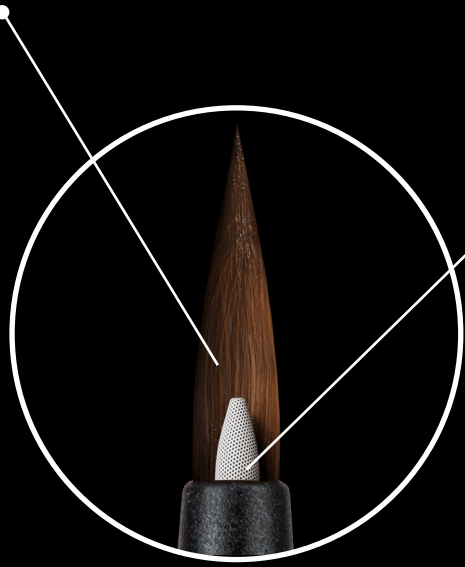

Absorbent Felt Nib



MPF Hybrid Brush Tips




Video

Patent pending (A.N. : CY202200003)

Due to the decline in the availability of high-quality Kolinsky hair, we have developed a solution that incorporates a felt nib insert within the hair bundle. This absorbent felt nib retains moisture, creating a consistently humid environment. As the hair bundle wraps around the felt nib, it helps maintain the hair in a continuously humid state.



Absorbent Felt Nib

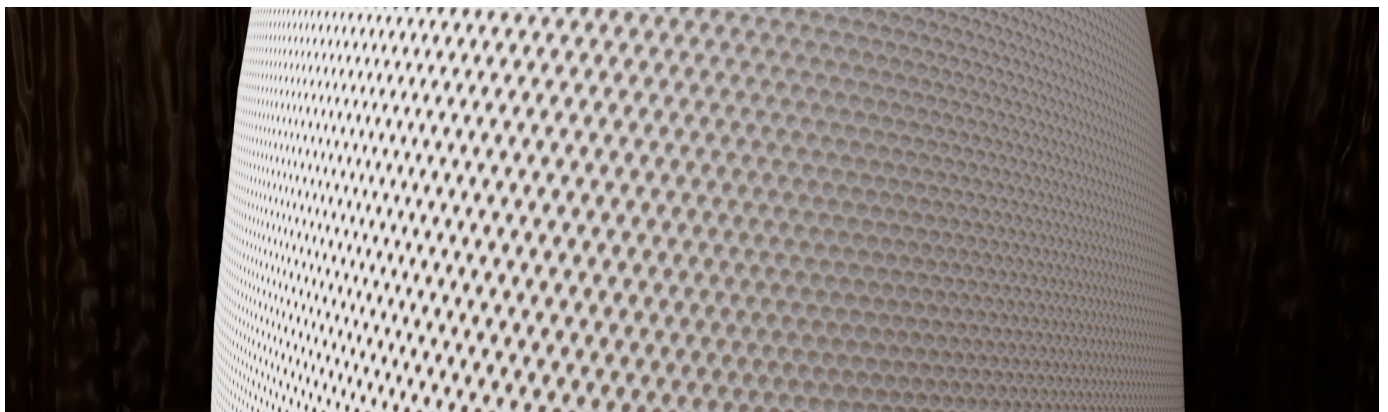

Compatible with every MPF Brush Handle.

70% Synthetic Hair

30% High grade Kolinsky

Absorbent Felt Nib



02

Composite Brushes & Tools




Video


MPF Composite Stain Brushes STRAIGHT

The MPF Composite Brush series features innovative straight brushes designed for precise application of composite stains. Crafted with a special Fiber optic filament, these brushes combine fine texture with firmness, making them ideal for intricate detailing.

The unique filaments provide a fine texture that allows for smooth and even application, ensuring flawless results.

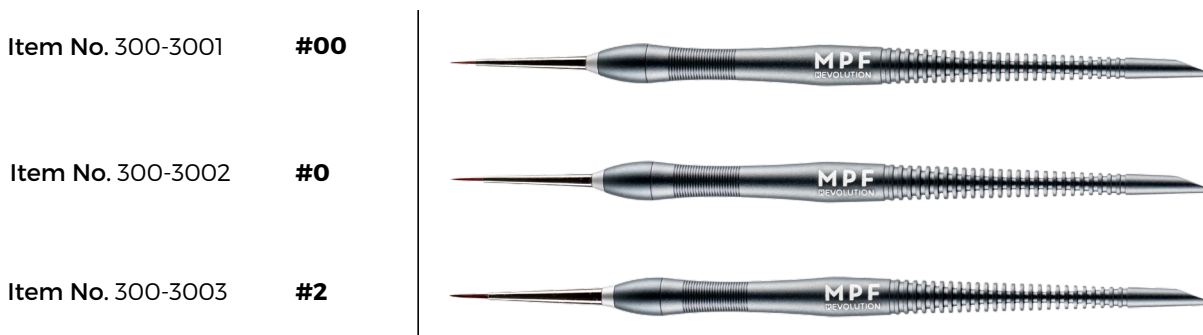

MPF Composite Stain Brushes ANGLED

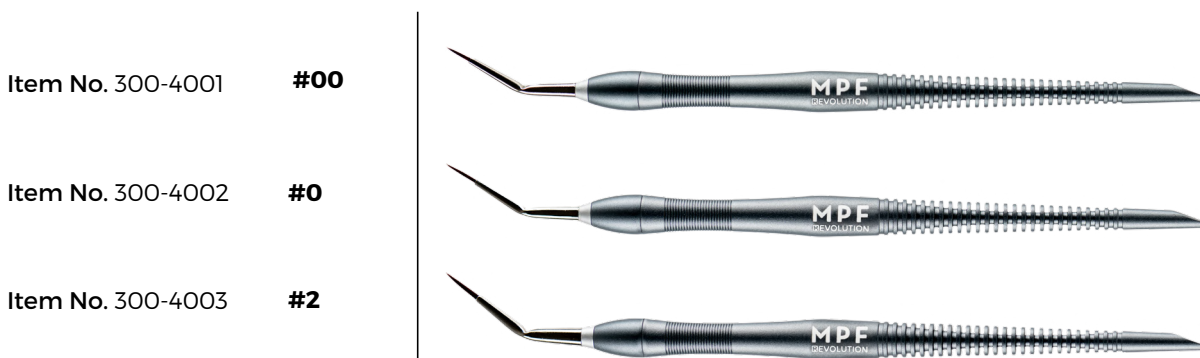

Item No. 300-3000

MPF Composite 3 in 1 Stain Brush Kit.
Straight



Item No. 300-4000


MPF Composite 3 in 1 Stain Brush Kit.
Angled

* For Replacement Tips, see page **26**

MPF Composite Stain (Straight) Replacement Tips



MPF Composite Stain (Angled) Replacement Tips

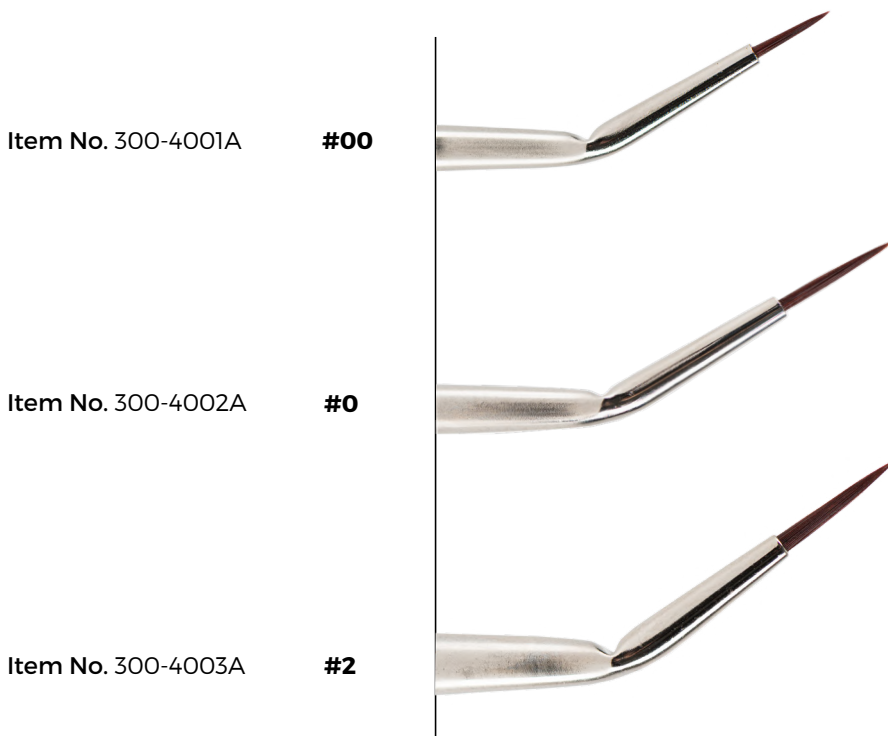



Video

MPF Composite Flat Brushes STRAIGHT

The MPF Composite Flat Brushes STRAIGHT are expertly designed for efficient application and detailing of composite materials. Their flat shape makes them ideal for broader strokes and precise coverage.

Made with durable, fine-textured Fibers that ensure consistent performance and superior results.

Item No. 300-5001 **#4**



Item No. 300-5002 **#6**



Item No. 300-5003 **#7**



MPF Composite Flat Brushes ANGLED

Item No. 300-6001 **#4**



Item No. 300-6002 **#6**



Item No. 300-6003 **#7**



Item No. 300-5000

MPF Composite 3 in 1 Flat Brush Kit.
Straight

Item No. 300-6000


MPF Composite 3 in 1 Flat Brush Kit.
Angled

* For Replacement Tips, see page **28**

MPF Composite Flat (Straight) Replacement Tips



MPF Composite Flat (Angled) Replacement Tips

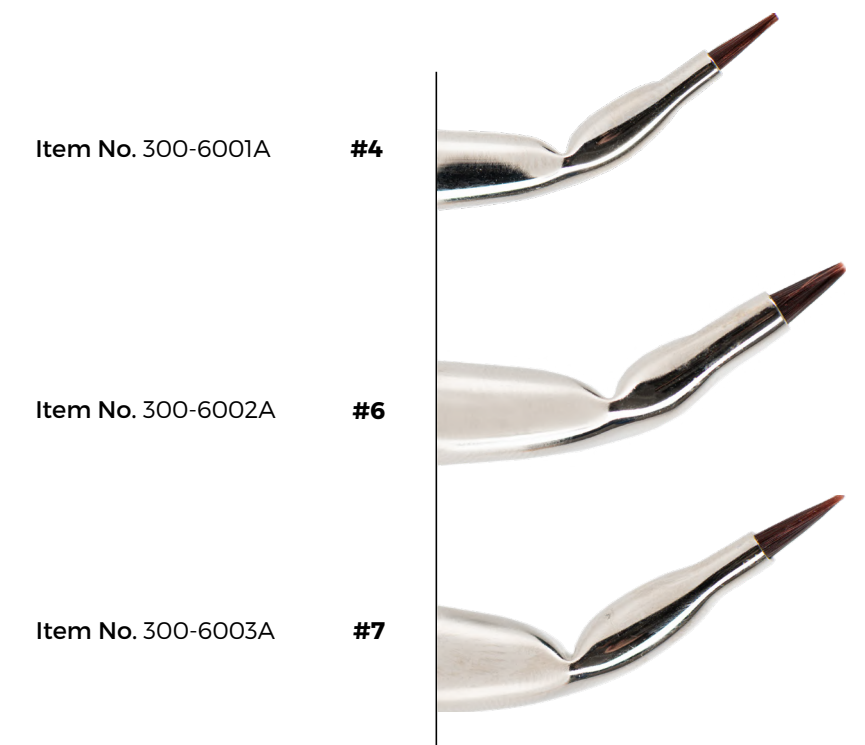

MPF Composite Application Tools (Direct)




Video

The MPF Composite Application Tools are expertly crafted to provide dental professionals with enhanced flexibility and precision for various restorative procedures. Designed for full mouth restorations, including gingival additions, single crowns, inlays, and onlays, these tools are essential for the application of dental composites.

Capable of creating stippling effects and intricate details, these tools are perfect for achieving high-quality restorations. Silicon tips added to the range for efficient shaping of larger volumes of material, ideal for posterior molars and gingival additions.

Custom create your Composite Application tool:

Step 1

Choose your handle  3 different colours for easy identification

Item No. 300-H-BL



Item No. 300-H-RG



Item No. 300-H-SL



Step 2

Choose your tips



#300-7001A

kolinksy stain brush



#300-7003B

non-stick silicone rounded point for gingival applications



#300-8001B

non-stick silicone angled tip for the shaping of larger areas of composite



#300-8003A

synthetic stippling brush for gingival characterizations



#300-7003A

large tapered spoon shape for general application of composite in all regions



#300-7002A

flat straight tip for anterior applications in the incisal regions



#300-8001A

medium sized synthetic flat brush for shaping and smoothening composite in the anterior region



#300-8002A

small tapered spoon for occlusal and incisal details



#300-7002B

angled tapered tip for interproximal applications



#300-8004B

angled spatula for mixing large qts of composite and for the shaping of gingival intended regions



#300-7001B

extra sharp deep fissure pointer for occlusal characterizations



#300-8004A

sharp deep fissure pointer for occlusal characterizations



#300-7004A

large angled occlusal groove creator



#300-7004B

long tapered deep groove and fissure creator



#300-8002B

angled deep fissure shaper for occlusal and anterior characterizations



#300-8003B

angled plunger for compressing large amounts of composite

MPF Composite Application Tools (Direct) Kits

Each Kit contains all 3 colors handles plus 6 specific tips.

Kit 01

#300-7000




Kit 02

#300-8000

Step 3

Assemble

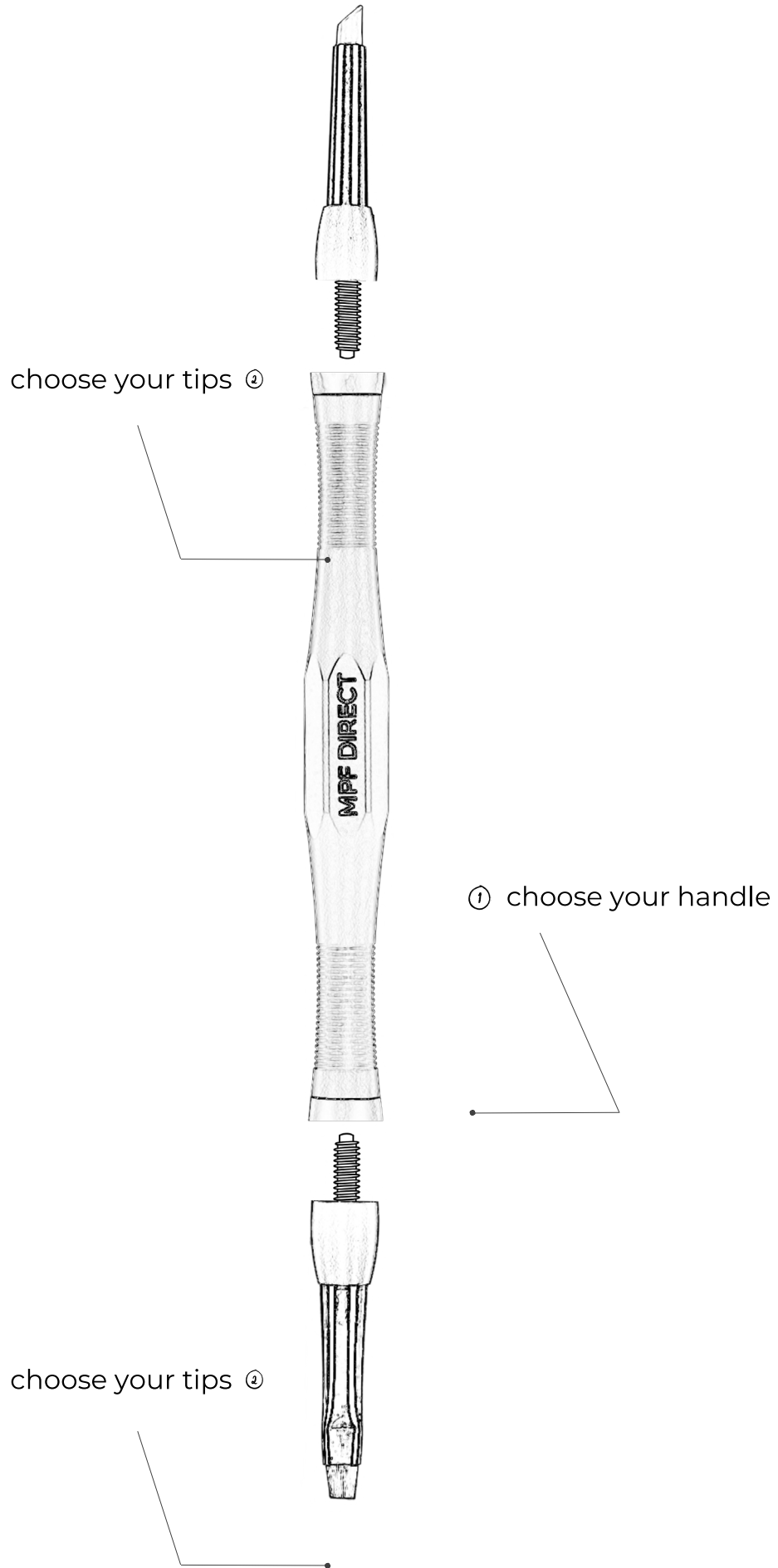

03

Palettes



Temperature Display Screen

Temperature Buttons



MPF
BRUSH C2
COMP TRAY

MPF COMPOSITE
WARMING PALETTE

TEMP

On/Off Switch

Power Inlet

Heating Tray

MPF Composite Warming Palette



Video

The MPF Composite Warming Palette is a groundbreaking tool designed to enhance the workability of composite materials, making them more malleable for precise application in dental restorations. Ideal for both small and large restorations, this palette is particularly beneficial in cases requiring extensive gingival adaptation. Made from Makrolon polycarbonate, the case provides stable UV protection to prevent premature polymerization of composites.


Item No. 115-9000

Includes: Universal power plug cable, blade scraper.

The lid can be positioned at any height, minimizing exposure to ambient light and maintaining optimal working conditions. The 6mm thick glass plate is heated from below, featuring an element connected to the main body, ensuring durability and resistance to breakage.

The palette features an On/Off button Located on the top right, it also allows you to switch the temperature display between Celsius and Fahrenheit. The LCD display shows the current temperature and allows adjustments using the + and - buttons. The palette requires a few minutes to reach the selected temperature.

The palette comes with a blade scraper for removing composite debris without scratching the glass surface. After scraping, the glass can be easily cleaned with a paper towel and alcohol.

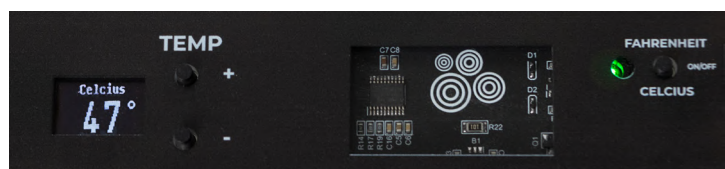

The palette will automatically turn off after 5 hours in case the unit has been left on by mistake.

Marble Ceramic Palette




Video

The Marble Ceramic Palette combines an Italian marble with our innovative design to provide an optimal workspace for ceramic mixing. This palette is specifically crafted to maintain the ideal consistency and workability of your ceramic materials.

Item No. 115-1000

Includes: marble tray, 2 worksheets, liquid dispenser bottle

The high-quality marble allows for controlled water permeation, keeping your ceramic mix consistently wet and workable. The sturdy construction ensures longevity and stability during use.

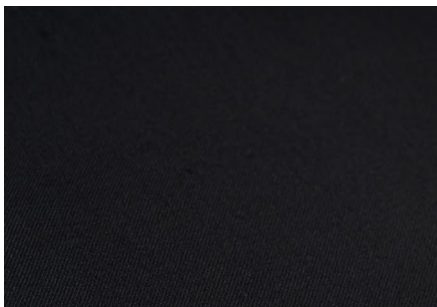

Item No. 115-1001

Nylon Palette Worksheets 2/pk



Item No. 115-1002

Marble Stone

Glass Ceramic Palette



Video


The Glass Ceramic Palette is designed for ceramists who prefer a smooth glass working surface. This palette utilizes an innovative system to ensure your ceramic mix remains consistently moist, enhancing your mixing and application processes.

Item No. 115-0000

Includes: glass tray, agate spatula, 4 strips, liquid dispenser bottle.

The Glass Work Surface provides a sleek and easy-to-clean area for mixing ceramics, allowing for precise control and application. The Base is built to withstand regular use while maintaining stability during work.

The Synthetic Leather Strips effectively keep the ceramic mix continuously wet. Water from the base permates through the strips, ensuring a homogeneous ceramic mix.

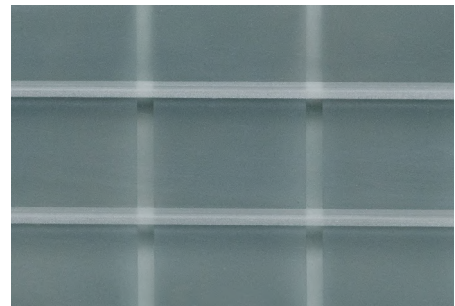

Item No. 115-0001

Palette Absorption Strips 4/pk

Item No. 115-0002

Replacement Glass



CAD Artistry Palette



Video


The CAD Artistry Palette is an innovative solution designed specifically for ceramists utilizing CAD-CAM technology. With dual working platforms, this palette provides optimal moisture and consistency control for ceramic materials.

Item No. 115-7000

Includes: tray, 2 strips, 1 worksheet, sponge, liquid dispenser bottle.

The CAD Artistry Palette was designed with two different options of humidification to keep the proper moisture and proper consistency of ceramic material. The two platforms on the **silicone tray** allow you to choose between humidification with water or with liquid. The new silicon is highly compatible with the modern staining glaze materials of today.

Choose between working on a nylon sheet or using the synthetic strips. Eleven stain compartments act as storage for your favourite stains and colours.

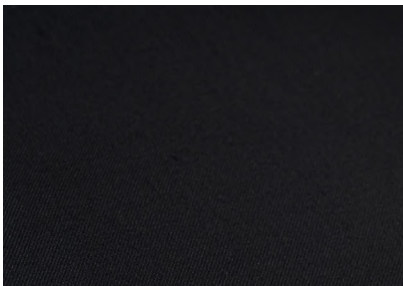


Item No. 115-7001

Nylon Palette Worksheet



Item No. 115-0001

Palette Absorption Strips 4/pk

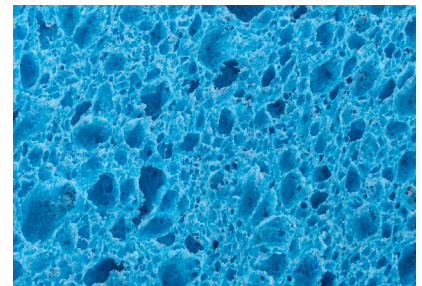

Item No. 115-4005

Replacement Sponge

Stain/Opaque Palettes



Video

The Stain/Opaque Palettes have been redesigned using silicone to better accommodate modern stain materials. This innovative palette offers smart storage solution for stains and opaque pastes, ensuring durability and convenience.

Item No. 115-5000
Black

Item No. 115-6000
White

The updated silicone construction is more compatible with contemporary stain materials, providing flexibility and resilience. Durable Plastic Structure built to withstand regular use, ensuring long-lasting performance. The tightly sealing lid prevents dust contamination and keeps materials safe from drying out, maintaining their quality over time.

The unique angled design allows for easy access to materials, enhancing usability and efficiency during application

A handy colour chart on the inside of the lid keeps your materials organized and easily accessible.



Evolution Master Ceramic Palettes



Video

The Evolution Master Ceramic Palette is designed for ceramists seeking a sophisticated and user-friendly workspace. With a sleek glass work surface and innovative features, these palettes ensure optimal mixing and application of ceramic materials

Item No. 115-4000

Black

Item No. 115-4100

White

Includes: glass tray, 4 strips, 1 nylon worksheet, sponge, stain brush

The glass work surface provides a smooth, durable area for mixing ceramics, allowing for precise control and easy cleanup.

The synthetic leather strips keep the ceramic mix continuously wet, ensuring optimal consistency and workability. The concealed stain palette drawer is a discreet storage solution for your favorite stains, keeping your workspace organized and clutter-free.



Item No. 115-4006

Replacement Glass Tray (Black)



Item No. 114-4007

Replacement Glass Tray (White)



Item No. 115-4001

Replacement Nylon Worksheet



Item No. 115-4003

Replacement Synthetic Strips 4/pk



Item No. 115-4005

Replacement Sponge



04

Paste Ceramic Characterization Kit

MPF Paste Ceramic Characterization Kit

MPF Paste Ceramic Characterization Kit is specifically designed for the latest thixotropic, ceramic-based materials. Providing superior control over flow and consistency. Each brush offers optimal stiffness and tip precision, allowing you to sculpt, layer, and detail paste ceramics with minimal effort and maximum accuracy.

Item No. 200-6000-M

Master Kit



Item No. 200-6000-M

Basic Kit



MPF Surface Striations Flat Brush Tip

Surface smoothing with controlled striation formation

Item No. 200-6001



MPF Structure Layering Brush Tip

Precision staining & micro-layering for paste ceramics

Item No. 200-6002



MPF Gingiva Structure Layering Brush Tip

Larger layering volume for gingival structure pastes

Item No. 200-6003



MPF Gingiva Silicone Shaper Tip

Gingiva characterization & morphology creation with paste ceramics

Item No. 200-6004



05

Green State Zirconia Brushes & Tools

MPF Zirconia Stain Brushes & Kit

Introducing the MPF Zirconia Stain Brushes, a revolutionary solution for all your staining needs. The high-quality synthetic bristles are specifically designed to work effectively with water-based liquids and acids which cause corrosion to metal and lacquer treated brushes. This series is made completely from recyclable plastic.

Available in a complete kit for all your staining requirements.

Item No. 105-1001

#1



Item No. 105-1002

#2



Item No. 105-1003

#3



Item No. 105-1000

MPF Zirconia Stain 3 in 1 Brush Kit.

MPF Zirco Texturizer



Video

The MPF Zirco Texturizer is specifically designed for creating soft, subtle textures on pre-sintered zirconia. This essential tool helps enhance the aesthetic quality of zirconia restorations before the sintering process.

Item No. 400-1000

Ideal for use before sintering, allowing you to achieve desired textural effects that remain after the process. Constructed from durable stainless steel, the brush effectively creates asymmetric lines that mimic natural texture. The texturizer also serves to remove any fine zirconia powder from the outer surface of the restoration, ensuring a clean finish.



MPF Zir Grip

The MPF Zir Grip Tweezers are expertly designed for working with all types of ceramic and zirconia restorations. These tweezers offer precision and functionality.



Video

Item No. 400-5000

Crafted for comfort and ease of use, allowing for extended handling without fatigue. The tips are made from sintered zirconia, specifically shaped to assist in the accurate positioning of crowns, veneers, and small bridges on refractory pins before firing. Unlike traditional metallic tweezers, the zirconia tips prevent the transfer of dark marks to the restoration, ensuring a clean and flawless finish.



MPF Ceramo Blade



Video

The MPF Ceramo Blade is an essential tool for creating intricate surface designs on pre-sintered zirconia. This precision instrument allows ceramists to achieve detailed effects that enhance the aesthetic quality of restorations.

Item No. 400-3000

Item No. 400-3001
Replacement Tip

Designed specifically for pre-sintered zirconia, the blade enables the creation of detailed features such as crack lines and deep fossa in the occlusal region. The blade's sharpness facilitates easy removal of zirconia material, allowing for the creation of wear facets and other creative details prior to final sintering. Ensures accurate and controlled carving, essential for high-quality restorations.



MPF Zirco Cleaner

The MPF Zirco Cleaner is specifically designed to simplify the cleaning process for zirconia milled crowns, bridges, and implant restorations. This tool is essential for maintaining the quality and aesthetics of zirconia work.




Video

Item No. 400-4000

Item No. 400-4001
Replacement Tip

The brush features tightly bundled synthetic hair, providing excellent stability and effectiveness in removing fine zirconia powder. The brush tip is tailored for navigating difficult-to-reach areas, ensuring thorough cleaning without leaving residues that could affect fitting. Designed to address the challenge of cleaning fine powder that is often hard to see, making it easier to achieve a pristine finish.

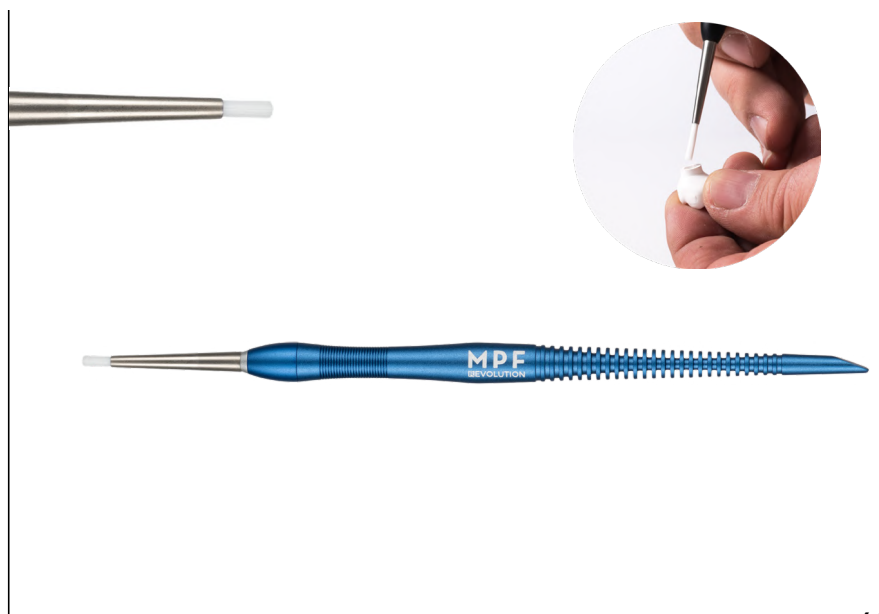

MPF Zirco Scorer




Video

The MPF Zirco Scorer is designed specifically for creating intricate surface textures on pre-sintered zirconia. This precision tool allows ceramists to achieve detailed effects that enhance the overall appearance and functionality of zirconia restorations.

Item No. 400-2000

Item No. 400-2001
Replacement Tip

The Zirco Scorer enables the crafting of intricate details, including pit fissures, grooves, and crack lines, providing a natural and realistic finish. Handling: Made from zirconia itself, the scorer prevents any contamination of the raw material, ensuring purity and quality in your work.

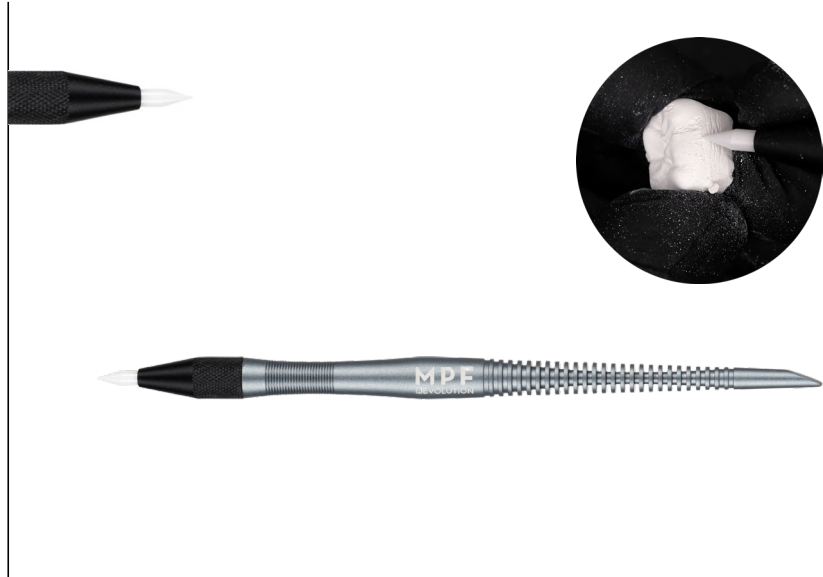

MPF Ceramomix™ Zirconia Spatula


The MPF Ceramomix™ Zirconia Spatula is a specialized ceramic mixing instrument designed to prevent contamination of ceramic materials while enhancing efficiency in your mixing process.



Video

Item No. 108-1000

The Spring Activated Technology design relieves pressure on the zirconia blade during mixing, significantly extending the lifespan of the spatula by tenfold. The spatula is engineered to ensure that ceramic materials remain uncontaminated, preserving the quality of your mixes. Built to withstand regular use, the spatula combines functionality with durability, making it a reliable tool in your toolkit.

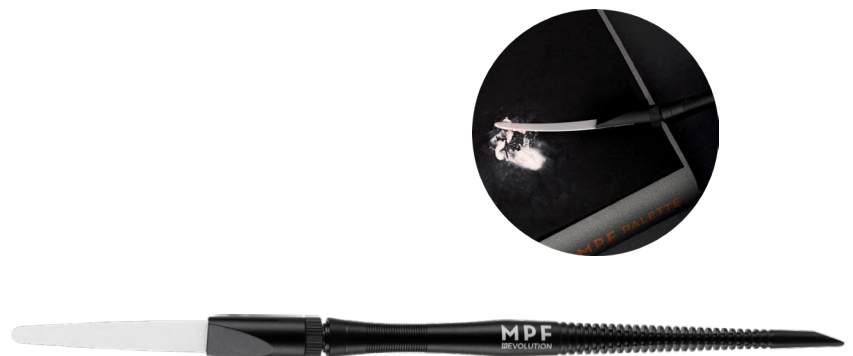

06

Lab Materials

MPF Veneer & Crown Stabilizing Paste

MPF Veneer & Crown Stabilizing Paste is a specialized adhesive designed to securely hold metal and metal-free prosthetics to model dies. Formulated to provide strong adhesion, ensuring that your prosthetics remain firmly in place during the fabrication process. The paste burns off cleanly without leaving any residue, making it ideal for use in dental applications where cleanliness is crucial.

Item No. 150-0001



MPF Peg Fix

MPF Peg Fix is an innovative, environmentally friendly refractory stabilizing paste designed for use with both metal and metal-free prosthetics. Made from PCW Fibers, this paste is a sustainable choice for dental professionals looking to reduce their environmental impact. Easily removed from your prosthetics with a short blast of compressed air, ensuring a clean surface without any residue left behind.

Item No. 150-0002
(1 tube)

Item No. 150-0003
(3 tubes)



MPF Polishing Kit

MPF Polishing Kit is designed as a complete, reliable solution for dental technicians, combining essential tools with carefully selected rotary accessories to streamline both pre-polishing and polishing steps. Each component is crafted to enhance efficiency and precision, allowing you to achieve smooth, high-gloss surfaces on every restoration.

Item No. 113-1000-K



MPF Texture Eyes



Item No. 153-0000

MPF Texture Eyes is a unique copper paste that aids in visualizing and replicating the morphology and structure of teeth.

MPF All-Bright Diamond Paste



Item No. 113-1000

MPF All-Bright Diamond Paste is a polishing solution that provides a lustrous, high-shine finish to ceramic and composite surfaces.

MPF Rainbow Drops Ceramic Colorants

MPF Rainbow Drops Ceramic Colorants are high-quality colorants specifically designed to enhance your ceramic mixes. These vibrant colorants allow for creative expression and customization in your ceramic artistry.

Item No. 151-0000 (Kit)

Item No. 151-0001 (Red)

Item No. 151-0002 (Blue)

Item No. 151-0003 (Orange)

Item No. 151-0004 (Green)

Item No. 151-0005 (Yellow)

Item No. 151-0006 (Violet)



MPF Debubblizer

The MPF Debubblizer is an advanced solution designed to enhance the quality of dental model casts by significantly reducing surface tension and eliminating bubble formation. Clinically designed to effectively minimize surface tension in impression materials, ensuring optimal results during casting. Simply mist the Debubblizer on the silicone, rubber base, or alginate impression materials to relieve surface tension, leading to flawless replications. Enhances the flow of gypsum and investment materials, ensuring they fill all impressions without obstruction.

Item No. 157-0000
300ml

Item No. 157-1000
1000ml


The MPF Debubblizer comes in a 300ml bottle including spray nozzle and also 1000ml refill bottles.



MPF Honeycombs & Pins

MPF Honeycombs & Pins are featuring innovative 'hot zone' technology, specifically designed for enhanced performance in crown and implant bound prosthesis baking within the hot zone of the furnace.

Item No. 114-1000
Honeycomb Tray with 10 Ceramic Pins

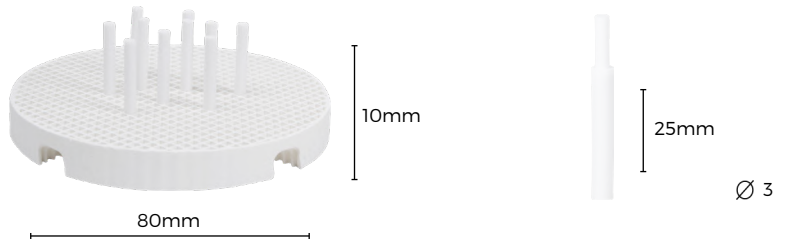

Item No. 114-1001
Ceramic Pins, 10 pcs

Item No. 114-2000
Honeycomb Tray 10 Metal Pins



Item No. 114-2001
Stainless Steel Pins, 10 pcs

Item No. 114-5000
Honeycomb Tray 10 Ceramic Pins




Item No. 114-5001
Ceramic Pins, 10 pcs (5 of each shape)

Item No. 114-3000
Honeycomb Tray for Implant bound CR & BR with 4 ceramic pins



Item No. 114-3001
Custom firing Ceramic Pins for Implant bound CR & CB, 4 pcs.

Item No. 114-4000
Honeycomb Tray for Veneers
40mm x 20mm

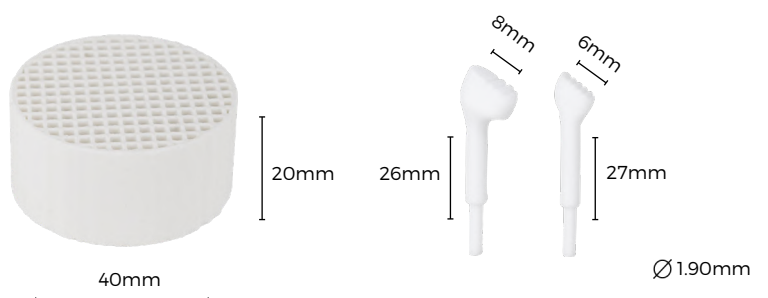

Item No. 114-4001
MPF Alumina Veneer Pins 6pcs

MPF UltraSep

MPF UltraSep has been specifically formulated to provide superior insulation for acrylic materials, ensuring optimal separation between acrylic components and 3D-printed models (FDM, FFF, DLP, SLP). Its advanced formulation creates an ultra-thin insulation layer with a single application, which can be effortlessly removed.

Compatible with most 3D printed resins. For external use on model surfaces only. Not a medical device.

Item No. 810-0000

1pcs

Item No. 810-0000K

3pcs

MPF UltraSep is ideal for dental laboratory procedures where insulation is required between methyl methacrylate-based acrylics and other surfaces, including:

- Separation of printed models during acrylic layering
- Insulation between acrylic modeling materials and lab-generated models
- Use in workflows involving custom trays, record bases, or provisional lab component



07

Lab Tools

MPF Brush Hydrator

The MPF Brush Hydrator is designed to keep your brushes perfectly hydrated throughout the day, ensuring optimal performance and longevity.

Item No. 109-0000

Simply pour water into the glass reservoir and mount your brush on the hydrator for about 5 minutes. This process allows the water to permeate the hair bundle.

ATTENTION: Hydrator is not a storage solution.



* brush not included

MPF Brush Pool

The MPF Brush Pool is an ergonomic and durable solution designed for the cleaning and storage of your favourite MPF brushes. The base is made from high-quality silicone, ensuring easy cleaning and fast drying after the use.

Item No. 111-0000
Brush Pool

Item No. 111-0001
Replacement Bowl

Item No. 111-0002
Replacement Sponge



MPF Master Klip

The MPF Master Klip Tweezers are an innovative, spring-activated tweezers designed for precision handling of delicate dental prosthesis. Ideal for holding crowns, veneers, and other delicate prosthetics during staining or while making additions, ensuring accuracy and control.

Item No. 112-1000

Made from durable stainless steel, the pointed holders ensure longevity and reliability while maintaining a firm grip on delicate items.



Platinum Foil Technique Tool Kit by Ricardo Schäfer

The Platinum Foil Technique Kit, has been developed in collaboration with renowned expert Ricardo Schäfer from Argentina, This kit is designed to streamline the process of working with platinum foil in dental applications.

Item No. 117-0000



MPF Drill Stands

The MPF Drill Stands are designed to provide efficient storage for your burs, keeping them organized and easily accessible based on your specific tasks. Equipped with non-slip silicone feet for enhanced grip, preventing movement during operation.



Item No. 155-0001



Item No. 155-0002



Item No. 155-0003



Item No. 155-0004

MPF Surgical Drill Stands

The MPF Surgical Drill Stand combines functionality and style, featuring a colourful anodized aluminium design and a three-level structure for optimal organization of your burs.

This innovative three-level design provides easy visibility of drill heads, with different levels accommodating various types of burs: **Lower Level:** Positions up to 26 friction fit burs, **Middle Level:** Holds up to 19 friction fit burs, **Upper Level:** Fits up to 12 regular hand motor burs.

The Surgical drill stand is connected to a flat base that allows the upper section to spin freely for quick and easy access to the appropriate bur for any task. Designed for infection control, the MPF Surgical Drill Stand is fully autoclavable, ensuring a sterile environment.

The Drill Stand should not be put into an ultrasonic with chemicals !

Optimum sterilisation process is 99% alcohol, using an ultrasonic and then an autoclave complete program.

Item No. 156-0000

Aqua




Item No. 156-0001

Citrus



Item No. 156-0002

Rosegold

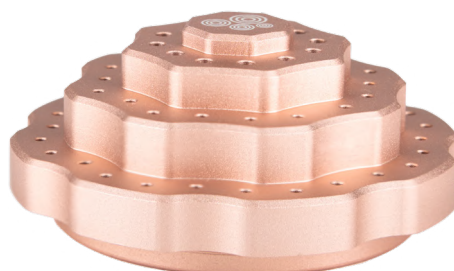

MPF Implant Lab Screw Drivers



Video

MPF Implant Drivers are ergonomically designed tools that provide precision and durability for dental implant procedures. The screwdriver bits are crafted from top-grade tool steel, heat-treated for maximum hardness and resilience. The handle is made from lightweight aluminium, providing comfort and ease of use during procedures.

Item No. 154-0001
Hexagon 0.9



Item No. 154-0002
Hexagon 0.95



Item No. 154-0003
Hexagon 1.2



Item No. 154-0004
Hexagon 1.27



Item No. 154-0005
Hexalobular T5.5



Item No. 154-0006
Hexalobular T6



Item No. 154-0007
Hexalobular Ball



Item No. 154-0008
Slot



Item No. 154-0000
Kit

Kit comes with an aluminum base for instant access and storage.
(see photo at page 63)



MPF IMPLANT DRIVERS

- 1. HEXAGON 0.9
- 2. HEXAGON 0.95
- 3. HEXAGON 1.2
- 4. HEXAGON 1.27
- 5. HEXALOBULAR T5.5
- 6. HEXALOBULAR T6
- 7. HEXALOBULAR BALL
- 8. SLOT



- 8. SLOT
- 7. HEXALOBULAR BALL
- 6. HEXALOBULAR T6
- 5. HEXALOBULAR T5.5
- 4. HEXAGON 1.27
- 3. HEXAGON 1.2
- 2. HEXAGON 0.95
- 1. HEXAGON 0.9



MPF Crown Holders



Video

MPF Crown Holders are specialized tools designed to securely hold crowns in place during dental procedures, ensuring safety and precision without damaging the workpiece. Made with a special gum material that protects crowns from chipping or marking while securely holding them in position.

Item No. 152-0001
3-4mm Handle

Item No. 152-0002
2-2.5mm Handle

Ensure precision and safety in your dental work with MPF Crown Holders, the ideal solution for handling crowns with care.



Item No. 152-0003
White

Item No. 152-0004
Black

The adhesive gum does not leave any sticky or greasy residue on crowns or fingers, maintaining cleanliness during procedures.



Item No. 152-0000
MPF Crown Holder
Complete Kit



MPF Spin Fix

The MPF Spin Fix is an innovative rotating table designed to enhance the workflow of dental ceramists. It features an alumina plate housed within the rotating upper plate, providing a stable and efficient platform for your crowns and bridges to cool down.

Item No. 106-0000

*honeycomb not included



MPF Bench Stand

The MPF Bench Stand is a versatile bench stand designed to keep your brushes and tools organized and within easy reach, enhancing your workflow during dental procedures.

Item No. 104-4000



Ceramic Table Mat

Ceramic Table Mats are an excellent addition to any dental laboratory, designed to enhance both aesthetics and functionality. These mats provide a beautiful backdrop for your workspace while offering practical benefits for technicians.



Item No. 110-0000



Item No. 110-0001

Liquid Dispensers



Item No. 115-0005
20 & 100 ml

Item No. 115-0003
Empty Modeling Liquid
Dropper 100ml

Item No. 115-0004
Empty Stain Liquid Dropper
20ml

MPF Cube Stand



Item No. 300-2000

Ergonomic stand for storing your brushes. Made from anodized aluminum.

Scribe™ Instrument by Joshua Polansky

The Scribe™ Instrument is a highly specialized tool designed for creating intricate surface textures and crack lines on baked ceramics and composite materials. Its sharp scoring capability allows for exceptional detail and precision, making it an invaluable addition to any ceramist's toolkit.

Item No. 107-0000

Item No. 107-0001
Replacement Tip

Ideal for use on both baked ceramics and composite materials, making it adaptable for various applications in your work.



MPF Agate Mixing Spatula

The MPF Agate Mixing Spatula is crafted from high-quality Brazilian agate stone, designed for mixing and extracting ceramic powders with ease and precision.

Item No. 108-0000

Made from Brazilian agate, known for its strength and safety, ensuring longevity and reliability in your mixing processes. The tail end is shaped at an angle, facilitating easy extraction of ceramic powder from bottles, reducing waste and enhancing efficiency. Specifically designed for all ceramic glaze pastes and other thixotropic materials, allowing for smooth and consistent results.



MPF Revolution Ceramic Sculpting Instruments

The MPF Revolution Ceramic Sculpting Instruments are crafted for precision and ease of use in ceramic carving in its wet stage. These tools are designed to meet the demands of ceramists seeking high-quality, reliable instruments for detailed work

Item No. 200-7001

Flat Blade



Item No. 200-7002

Angled Blade



Item No. 200-7003

Dagger



Item No. 200-7004

Lingual Spoon




Item No. 200-1701

Rplc Tip




Item No. 200-1702

Rplc Tip

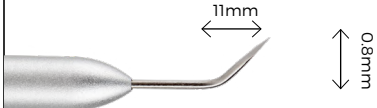


Item No. 200-1703

Rplc Tip

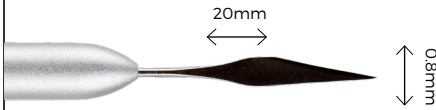

Item No. 200-1704

Rplc Tip

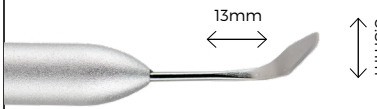

Item No. 200-7000

MPF Revolution Sculpting Instrument Kit

MPF Waxer Instruments

The MPF Waxer Instruments are versatile tools designed for dental professionals, featuring lightweight aluminium handles and replaceable stainless-steel blades for optimal performance.

Item No. 119-0001
Spoon/Fine Dropper.
Red



Item No. 119-0002
Large LeCron/Angled Spade.
Blue



Item No. 119-0003
Small LeCron/Large Dropper.
Purple



Item No. 119-1001
Spoon



Item No. 119-1002
Fine Dropper



Item No. 119-1003
Large LeCron



Item No. 119-1004
Angled Spade



Item No. 119-1005
Small LeCron



Item No. 119-1006
Large Dropper



Item No. 119-2000
3 in 1 Kit



08

Rotary Instruments

MPF ZMAX



Video

ZMAX diamond-impregnated stones and discs are expertly designed for the abrasion of sintered zirconia, delivering exceptional performance with a focus on heat management. Ideal for shaping and finishing sintered zirconia, these tools provide a smooth and efficient cutting experience.

RPM 10.000 - 14.000

Item No. 120-0001

Large Disc 22 x 4.5 mm



Item No. 120-0002

Small Disc 15 x 3.5 mm



Item No. 120-0003

Disc 22 x 2 mm



Item No. 120-0004

Thin Disc 22 x 0.7 mm



Item No. 120-0005

Small Barrel 4.8 x 13 mm



Item No. 120-0006

Large Barrel 6.5 x 13 mm



Item No. 120-0007

Cone 3.5 x 11 mm



Item No. 120-0008

Inverted Cone 12 x 6 mm



Item No. 120-0009

Inverted Small Cone 6 x 8 mm



MPF Zirco Cera GSS



The MPF Zirco Cera GSS Series rotary instruments are specifically designed for optimal performance in working with zirconia and ceramic materials. They provide exceptional cutting abilities while minimizing the risk of damage to the substructure.

RPM 10.000 - 14.000

Item No. 121-0001

Large Disc 22 x 4.5 mm



Item No. 121-0002

Small Disc 15 x 3.5 mm




Item No. 121-0003

Disc 22 x 2 mm



Item No. 121-0004

Thin Disc 22 x 0.7 mm

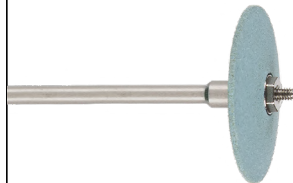

Item No. 121-0005

Small Barrel 4.8 x 13 mm



Item No. 121-0006

Large Barrel 6.5 x 13 mm



Item No. 121-0007

Cone 3.5 x 11 mm



Item No. 121-0008

Inverted Cone 12 x 6 mm



Item No. 121-0009

Inverted Small Cone 6 x 8 mm





MPF White Pearl Series

The MPF White Pearl Series consists of high-performance abrasive stones designed to efficiently tackle a wide range of materials, including ceramics, zirconia, and lithium disilicate.

RPM 10.000 - 14.000

Item No. 140-0001
Large Disc 22 x 4.5 mm



Item No. 140-0002
Small Disc 15 x 3.5 mm



Item No. 140-0003
Disc 22 x 2 mm



Item No. 140-0004
Thin Disc 22 x 0.7 mm



Item No. 140-0005
Small Barrel 4.8 x 13 mm



Item No. 140-0006
Large Barrel 6.5 x 13 mm



Item No. 140-0007
Cone 3.5 x 11 mm



Item No. 140-0008
Inverted Cone 12 x 6 mm



Item No. 140-0009
Inverted Small Cone 6 x 8 mm



MPF Rainbow Ceramic Polishers



Video

The MPF Rainbow Cutters are premium silicon and diamond-impregnated abrasive wheels designed specifically for use on ceramic and zirconia restorations. To achieve a lustrous finish on ceramic or zirconia, it's recommended to utilize all grades of the Rainbow collection in sequence.

Item No. 122-0001**#Coarse (Red)**

RPM 8.000 - 12.000

For minimal reduction on zirconia marginal regions and all types of ceramic materials.

**Item No. 122-0002****#Medium (Blue)**

RPM 10.000-12.000

For minimal reduction on zirconia marginal regions and for contouring all types of ceramic materials.

**Item No. 122-0003****#Fine (Grey)**

RPM 12.000-15.000

For the removal of minor scratches on all types of ceramic materials.

**Item No. 122-0004****#High Shine (Yellow)**

RPM 15.000 - 20.000

For the final high shine procedure on all types of ceramic and zirconia materials.




MPF Aluminium Oxide Abrasive Discs

MPF Aluminium Oxide Abrasive Discs are designed for durability and efficiency, specifically tailored for alloy restoration applications. Constructed from high-quality aluminium oxide, these discs are built to withstand the rigors of alloy restoration tasks.

Item No. 123-0002

Large (38 x 0.6 x 1.8) 100 pcs

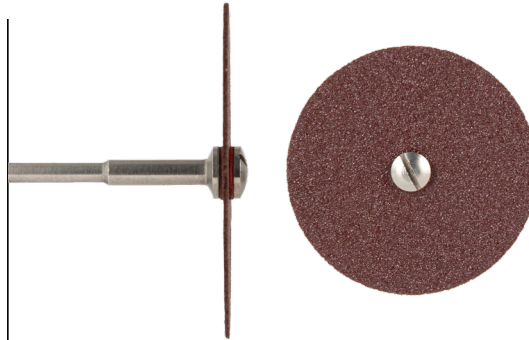


Item No. 123-0001

Middle (25 x 0.6 x 1.8) 100 pcs



Item No. 123-0003

Small (22 x 0.2 x 1.8) 100 pcs

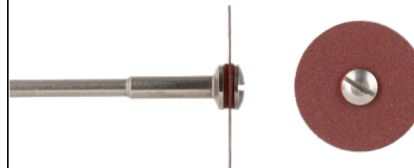

MPF Polishing Brushes & Polishing Mops

Item No. 124-0001 x100pcs

Item No. 124-0002 x50pcs

Item No. 124-0003 x25pcs



Item No. 125-0001 x100pcs

Item No. 125-0002 x50pcs

Item No. 125-0003 x25pcs



MPF Diamond Discs. **Mounted**

The MPF Diamond Discs have been developed especially for the use on all alloys, composite materials, ceramic materials and zirconia. The galvanically sintered edge offers maximum stability and a high diamond density. The serrated disc creates a cooling effect and they come in a medium and fine grit.

Replacements (5pcs) do not include the mount.

Item No. 141-0000 (1 pcs)

Item No. 141-0001 (5 pcs RPLC)

Fine (19mm x 0.13mm)

Item No. 142-0000 (1 pcs)

Item No. 142-0001 (5 pcs RPLC)

Medium (19mm x 0.13mm)

Item No. 143-0000 (1 pc)

Item No. 143-0001 (5 pcs RPLC)

Fine (22mm x 0.15mm)

Item No. 144-0000 (1 pcs)

Item No. 144-0001 (5 pcs RPLC)

Medium (22mm x 0.15mm)

Item No. 145-0000 (1 pcs)

Item No. 145-0001 (5 pcs RPLC)

(22mm x 0.20mm)




MPF Diamond Discs. **Unmounted**


The 22mm, 35mm & 45mm discs are perfect for the preparation of Gypsum dies, night guard materials and orthodontic aligners.



Item No. 146-0001
(22mm x 0.20mm) 1pc

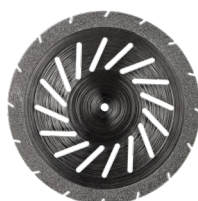

Item No. 147-0001
(35mm x 0.30mm) 1pc

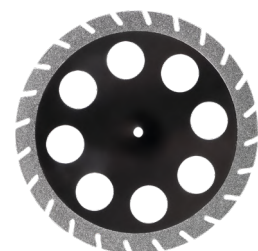

Item No. 148-0001
(45mm x 0.30mm) 1pc

Rubber Wheels & Points. **Mounted K Series.**

Expertly designed for cutting and polishing acrylic resin base plates and denture bases, making them essential tools for dental laboratories and clinical settings.

Item No. 133-0003 (x6pcs)



Item No. 133-0001 (x25pcs)

Mounted Rubber, coarse, diam/height: 13/20

Item No. 133-0012 (x6pcs)



Item No. 133-0010 (x25pcs)

Mounted Rubber, coarse, diam/height 10/24

Item No. 133-0021 (x6pcs)



Item No. 133-0019 (x25pcs)

Mounted Rubber, coarse, diam/height 15/18

Item No. 133-0006 (x6pcs)



Item No. 133-0004 (x25pcs)

Mounted Rubber, medium, diam/height: 13/20

Item No. 133-0015 (x6pcs)



Item No. 133-0013 (x25pcs)

Mounted Rubber, medium, diam/height: 10/24

Item No. 133-0024 (x6pcs)



Item No. 133-0022 (x25pcs)

Mounted Rubber, medium, diam/height: 15/18

Item No. 133-0009 (x6pcs)



Item No. 133-0007 (x25pcs)

Mounted Rubber, fine, diam/height: 13/20

Item No. 133-0018 (x6pcs)



Item No. 133-0016 (x25pcs)

Mounted Rubber, fine, diam/height: 10/24

Item No. 133-0027 (x6pcs)



Item No. 133-0025 (x25pcs)

Mounted Rubber, fine, diam/heightL 15/18



x6pcs
x25pcs
R.P.M 8000



x6pcs
x25pcs
R.P.M 10000



x6pcs
x25pcs
R.P.M 8000



x6pcs
x25pcs
R.P.M 8000



x6pcs
x25pcs
R.P.M 10000



x6pcs
x25pcs
R.P.M 8000



x6pcs
x25pcs
R.P.M 8000



x6pcs
x25pcs
R.P.M 10000



x6pcs
x25pcs
R.P.M 8000

Rubber Wheels & Points. **Unmounted F Series.**

The Unmounted F Series Rubber Wheels and Points are designed for precision cutting and polishing of semi-precious and precious alloys, making them essential tools for dental laboratories.

Item No. 130-0004 A (x100pcs)

Item No. 130-0004 B (x30pcs*)

Rubber Wheel, coarse, diam/height: 22/3



Item No. 130-0007 A (x100pcs)

Item No. 130-0007 B (x30pcs*)

Rubber Wheel, coarse, diam/height: 22/4



Item No. 130-0001 A (x100pcs)

Item No. 130-0001 B (x30pcs*)

Rubber Wheel, coarse, diam/height: 6.5/24



Item No. 130-0005 A (x100pcs)

Item No. 130-0005 B (x30pcs*)

Rubber Wheel, medium, diam/height: 22/3



Item No. 130-0008 A (x100pcs)

Item No. 130-0008 B (x30pcs*)

Rubber Wheel, medium, diam/height: 22/4



Item No. 130-0002 A (x100pcs)

Item No. 130-0002 B (x30pcs*)

Rubber Wheel, medium, diam/height: 6.5/24



Item No. 130-0006 A (x100pcs)

Item No. 130-0006 B (x30pcs*)

Rubber Wheel, fine, diam/height: 22/3



Item No. 130-0009 A (x100pcs)

Item No. 130-0009 B (x30pcs*)


Rubber Wheel, fine, diam/height: 22/4



Item No. 130-0003 A (x100pcs)

Item No. 130-0003 B (x30pcs*)

Rubber Wheel, fine, diam/height: 6.5/24

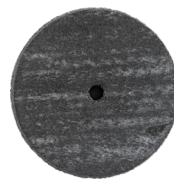


A: x100pcs
B: x30pcs
R.P.M 10000



A: x100pcs
B: x30pcs
R.P.M 10000



A: x100pcs
B: x30pcs
R.P.M 10000

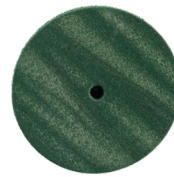

A: x100pcs
B: x30pcs
R.P.M 10000




A: x100pcs
B: x30pcs
R.P.M 10000



A: x100pcs
B: x30pcs
R.P.M 10000



A: x100pcs
B: x30pcs
R.P.M 10000

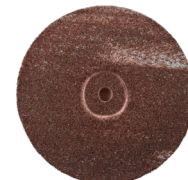

A: x100pcs
B: x30pcs
R.P.M 10000



A: x100pcs
B: x30pcs
R.P.M 10000

* 25pcs + 5 for free

* For Mandrels, see page **82**

Rubber Wheels & Points. **Unmounted H Series.**

The Unmounted H Series Rubber Wheels and Points are specifically designed for cutting and polishing dental ceramics, making them essential tools in dental laboratories.

Item No. 132-0004 A (x100pcs)

Item No. 132-0004 B (x30pcs*)

Rubber Wheel, coarse, diam/height: 22/3



Item No. 132-0007 A (x100pcs)

Item No. 132-0007 B (x30pcs*)

Rubber Wheel, coarse, diam/height: 22/4



Item No. 132-0001 A (x100pcs)

Item No. 132-0001 B (x30pcs*)

Rubber Wheel, coarse, diam/height: 6/22



Item No. 132-0005 A (x100pcs)

Item No. 132-0005 B (x30pcs*)

Rubber Wheel, medium, diam/height: 22/3



Item No. 132-0008 A (x100pcs)

Item No. 132-0008 B (x30pcs*)

Rubber Wheel, medium, diam/height: 22/4



Item No. 132-0002 A (x100pcs)

Item No. 132-0002 B (x30pcs*)

Rubber Wheel, medium, diam/height: 6/22



Item No. 132-0006 A (x100pcs)

Item No. 132-0006 B (x30pcs*)

Rubber Wheel, fine, diam/height: 22/3



Item No. 132-0009 A (x100pcs)

Item No. 132-0009 B (x30pcs*)


Rubber Wheel, fine, diam/height: 22/4



Item No. 132-0003 A (x100pcs)

Item No. 132-0003 B (x30pcs*)

Rubber Wheel, fine, diam/height: 6/22

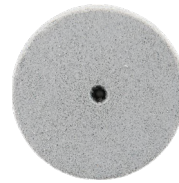


A: x100pcs
B: x30pcs
R.P.M 10000



A: x100pcs
B: x30pcs
R.P.M 10000



A: x100pcs
B: x30pcs
R.P.M 20000

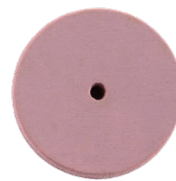


A: x100pcs
B: x30pcs
R.P.M 10000



A: x100pcs
B: x30pcs
R.P.M 10000



A: x100pcs
B: x30pcs
R.P.M 20000

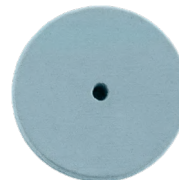

A: x100pcs
B: x30pcs
R.P.M 10000



A: x100pcs
B: x30pcs
R.P.M 10000




A: x100pcs
B: x30pcs
R.P.M 20000

* 25pcs + 5 for free


* For Mandrels, see page **82**

Tungsten Carbide Burs

Coarse

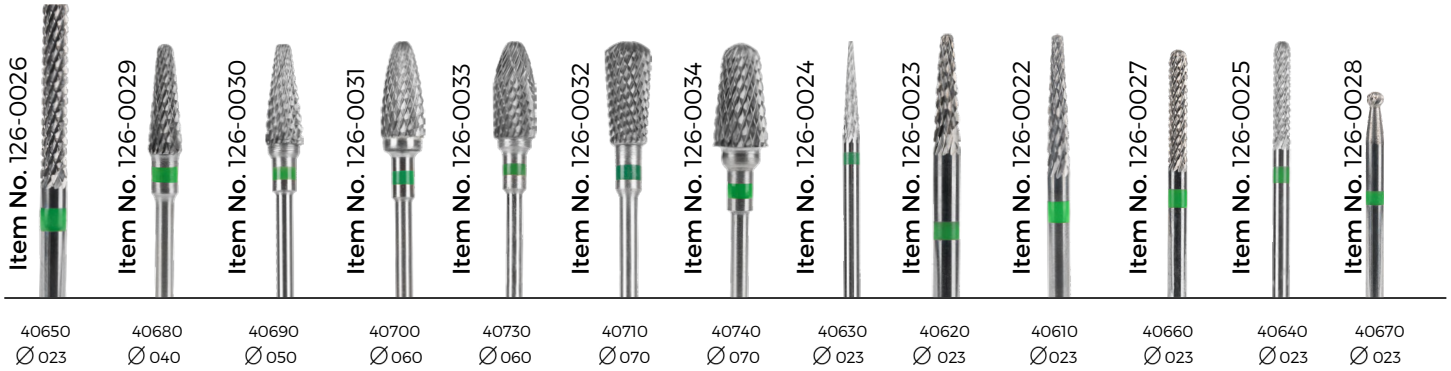


Super Coarse

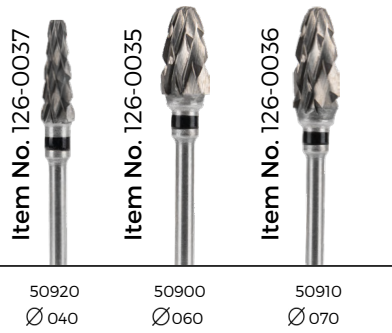

Fine



Super Fine

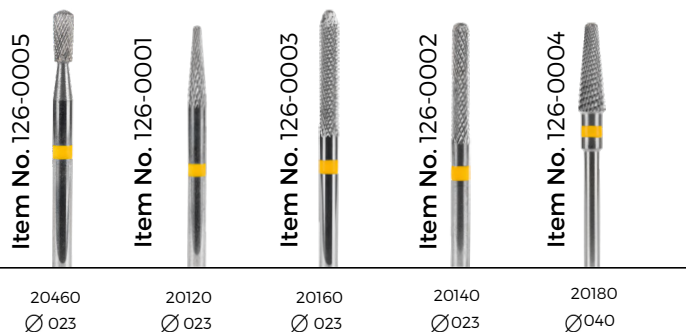


*HEAD DIAMETER: Ø

F.G Diamond Burs



Medium to Coarse

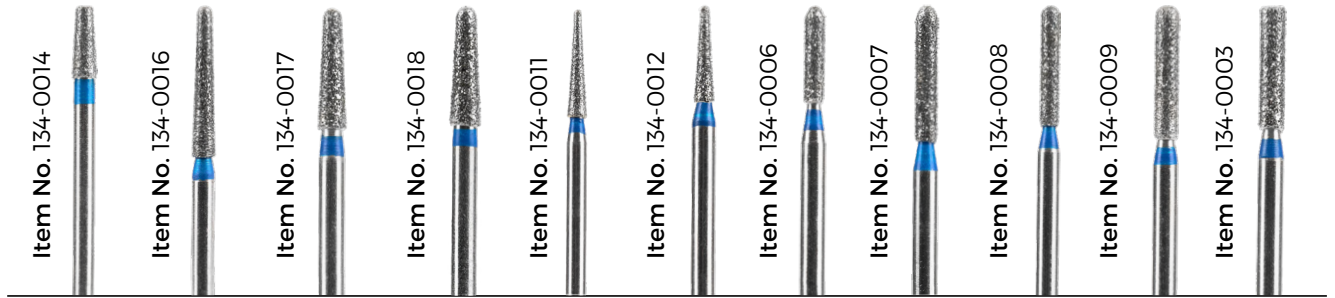

ISO No.	170/018	199/016	197/016	197/018	160/014	160/012	200/012	289/014	140/012	141/014	110/014
Product ID	TF-23	TR-12	TR-21	TR-24	TC-21	TC-26	SO-20	SO-21	SR-11	SR-12	SF-21

Medium to Coarse



ISO No.	257/018	111/012	111/014	234/014	237/010	237/018	010/012	010/018	001/012	001/016	068/042
Product ID	FO-32	SF-11	SF-12	EX-31	EX-41	EX-20	SI-46	SI-49	BR-46	BR-40	WR-13W

Medium to Fine



ISO No.	173/016	199/016	199/018	160/016	298/014	298/016	237/021	196/020	197/012	254/017	108/008	138-007
Product ID	TF-12F	TR-11F	TR-26F	TC-11F	FO-21F	FO-22F	EX-21F	CR-11F	CR-21F	FO-30F	CD-58F	CD-52F

* For Converters, see page **82**

MPF Corkscrew Mandrel

Stainless Steel high speed corkscrew mandrel for felt points & felt cylinders.

Item No. 128-0006 (6pcs)

Item No. 128-0015 (15pcs)

Shank size 2.35mm



Ø 5.0mm

MPF HP Mandrel

Stainless Steel high speed mandrel

Item No. 127-0006 (6pcs)

Item No. 127-0015 (15pcs)

Shank size 2.35mm



Ø 5.0mm

MPF HP to FG Converter

Stainless Steel mandrel converter from HP into FG

Item No. 135-0001 (1pc)

Item No. 135-0005 (5pcs)

Shank size 2.35mm



Ø 5.0mm

09

Dental Training School Kit

MPF Dental Training School

The MPF Dental Training School Kit is designed specifically for dental clinical students and dental technology students, providing them with essential tools to support their learning and future careers in dentistry.

Item No. 160-0001

MPF Lessman Knife Large

Large & Small wax up knife for all waxing applications, with a knife shape and also a tapered scoop shape in the one handle.



Item No. 160-0002

MPF Lessman Knife Small



Item No. 160-0003

MPF G.S Wax Knife

Plaster knife with a large anterior blade for cutting and a leverage shape tip at the other end of the handle, usually used when deflasking dentures from muffles.



Item No. 160-0004

MPF Plaster Spatula

Straight edge plaster spatula for all mixing purposes.



Item No. 160-0005

MPF Cement Spatula Large

Large & Small cement and acrylic mixing spatulas with a tapered tip and a flatter shaped rear tip, the handles have excellent grip due to the knurled surface of the handles.



Item No. 160-0006

MPF Cement Spatula Small



MPF Dental Training School

The kit includes 27 high-quality hand tools that are essential for various dental procedures and techniques.

Item No. 160-0007

MPF Mini Blade Holder

Microblade holder with replacement blades used for removing bubbles from plaster models and general dental model preparation.



Item No. 160-0008

MPF Tweezer with Lock

Contra angle locking tweezers that can be used for holding sutures, carbon paper for bite registration purposes, general use in a Clinical or laboratory situation.



Item No. 160-0009

MPF Scapel Handle

Number 3 Scapel handle that is used with several blade sizes for a variety of applications in a laboratory or clinical situation.



Item No. 160-0016

MPF Scissors

Straight scissors for the cutting of sutures, platinum foil and an abundance of applications in a clinical and laboratory situation.



Item No. 160-0017

MPF Decimal Caliper

Measuring callipers used for measuring the thickness of ceramic veneers and or metal framework, applied in both a clinical and laboratory situation.



MPF Dental Training School

Item No. 160-0010**MPF Concave Plier**

Concave orthodontic pliers that can be used for several applications such as creating several types of Orthodontic clasps in a Laboratory or Clinical situation.

**Item No. 160-0011****MPF Adam Plier**


Adam pliers orthodontic pliers that can be used for several applications such as creating several types of Orthodontic clasps in a Laboratory or Clinical situation.

**Item No. 160-0014****MPF Needle holder Straight**

Needle holder Straight can be used in surgical procedures , but also used in a dental laboratory for several applications.

**Item No. 160-0015****MPF Needle holder Curved**

Needle holder Curved can be used in surgical procedures , but also used in a dental laboratory for several applications.

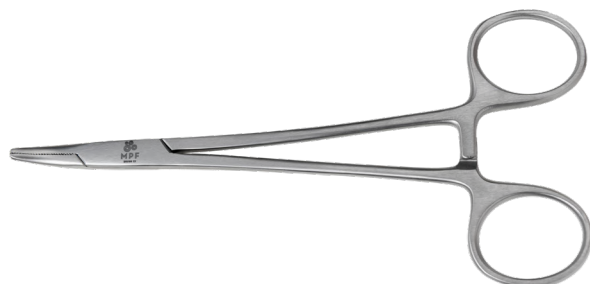

MPF Dental Training School

Item No. 160-0018

Alcohol Burner

Alcohol burner is used primarily in a clinical situation when working with wax or compound materials in a surgical situation.



Item No. 160-0019

Alcohol Torch

Alcohol torch can be used for the smoothing of wax surfaces in a clinical or laboratory situation .



Item No. 160-0020

Plaster Bowl

Silicon Plaster or Alginate mixing bowl which can be applied in both a clinical and laboratory situation.



Item No. 160-0021

Silicone Bowls

Silicon Bowls in small, medium and large sizes are primarily used for the mixing of powder to liquid materials which are applied in both a clinical and laboratory situation.



L



M



S



MPF
MARTIN

MPF Dental Training School Kit

Item No. 160-0000

Dental Training School Kit
*contains **27** tools in total.

Also included:

Combination Drill Kit

Mixing Plate.

MPF Classic Kolinsky Brushes #4 #6
(105-0004/105-0006).

Angle Spade/Large Lecron
(119-0002).

Large Dropper/ Small Lecron
(119-0003).

2 MPF Application Tools (Direct).



Premium Dental Instruments

MPF BRUSH Co.



MPF Brush Co.

web: www.mpbrush.com

tel: +357 24623314
